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NeoSoft, LLC

NEOSOFT

NS-03-039-0014 Rev 2
Copyright 2020 NeoSoft, LLC
All rights reserved

suiteHEART® Software

Addendum to Instructions for Use

NeoSoft, LLC

The logo for NeoSoft, featuring the word "NEO" in a bold, cyan-colored sans-serif font, followed by "SOFT" in a grey-colored sans-serif font. The letters are closely spaced and have a slight shadow effect.

NS-03-039-0014 Rev 2
Copyright 2020 NeoSoft, LLC
All rights reserved

Revision History

Rev	Date	Description of Change	Safety Related Update (Yes/No)
1	26APR2019	Updated for 5.0.0 product release. English manual previously NS-03-015-0006, Rev 7. Added French, German, Greek, and Italian languages which were previously: NS-03-015-0007 FR Rev 7 NS-03-015-0008 DE Rev 7 NS-03-038-0009 EL Rev 1 NS-03-038-0018 IT Rev 1	No
2	12NOV2019	Updated for 5.0.1 product release.	No

Manufacturer

NeoSoft, LLC
N27 W23910A Paul Road
Pewaukee, WI 53072 USA

Phone: 262-522-6120
website: www.neosoftllc.com

Sales: orders@neosoftmedical.com
Service: service@neosoftmedical.com

The Instructions for Use for this device are provided electronically in Portable Document Format, (.pdf). A pdf viewer is required to view the Instructions for Use. A paper copy of the Instructions for Use can be provided upon request, at no cost, within 7 calendar days, by emailing service@neosoftmedical.com.

Instructions for Use may be accessed in the following ways:


1. After launching the application, click "Help" or "About" from the main screen. Select the "Instructions for Use" option. The Instructions for Use will open in a pdf viewer.
2. If the original installation package received from NeoSoft is available, open the zip file and navigate to the "Documentation" folder, then the "Instructions for Use" folder and double click on the Instructions for Use.pdf in your language, denoted by EN - English, FR - French, DE - German, EL - Greek, IT - Italian, LT - Lithuanian, ES - Spanish, SV - Swedish, TR - Turkish, RO - Romanian, NL - Dutch.
3. Navigate to the folder where the application is installed. Locate the "Documentation" folder, open the "Instructions for Use" folder and double click on the Instructions for Use.pdf in your language, denoted by EN - English, FR - French, DE - German, EL - Greek, IT - Italian, LT - Lithuanian, ES - Spanish, SV - Swedish, TR - Turkish, RO - Romanian, NL - Dutch.
4. Electronic copies of Instructions for Use are also available at www.neosoftllc.com/neosoft/product_manuals/.

Medical device directive

This product conforms with the requirements of council directive 93/42/EEC concerning medical devices when it bears the following CE mark of Conformity:

European Representatives:

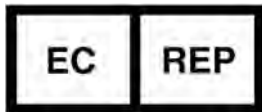EMERGO EUROPE
Prinsessegracht 20
2514 AP The Hague
The Netherlands

Canada:

Health Canada device license number: 99490


CAUTION: *United States federal law restricts this device to sale by, or on the order of, a physician.*

System Requirements

*Recommended video card(s); NVIDIA Quadro P4000 or P5000

suiteHEART® Software Minimum System Requirements (Multiuser Environment)
Supported Operating Systems: <ul style="list-style-type: none"> • Windows 2012 R2 Standard, Windows Server 2016 Standard, Edition (64 bit versions) or Windows Server 2019 Standard
Citrix XenApp 7.7 or higher
Microsoft® .NET Framework 4.0
Video card*(NVIDIA); 8GB RAM (minimum), minimum of 1664 CUDA Cores, support for OpenGL 4.0 and support for 1920 x 1080 resolution or higher
Monitor shall support 1920 x 1080 resolution or higher
Anti-Virus/Malware software installed
Minimum 64GB memory
Available hard drive space minimum 500GB
PDF viewer, Adobe Reader 11.0 or higher
An open network port for DICOM transfer over network
Static IP address (recommended)
Solid-State hard drive(s)
Intel® Core™ processor(s) (4 cores)
Citrix XenApp client displays must not exceed scaling/zoom levels greater than 125%

Citrix Installation (Example)

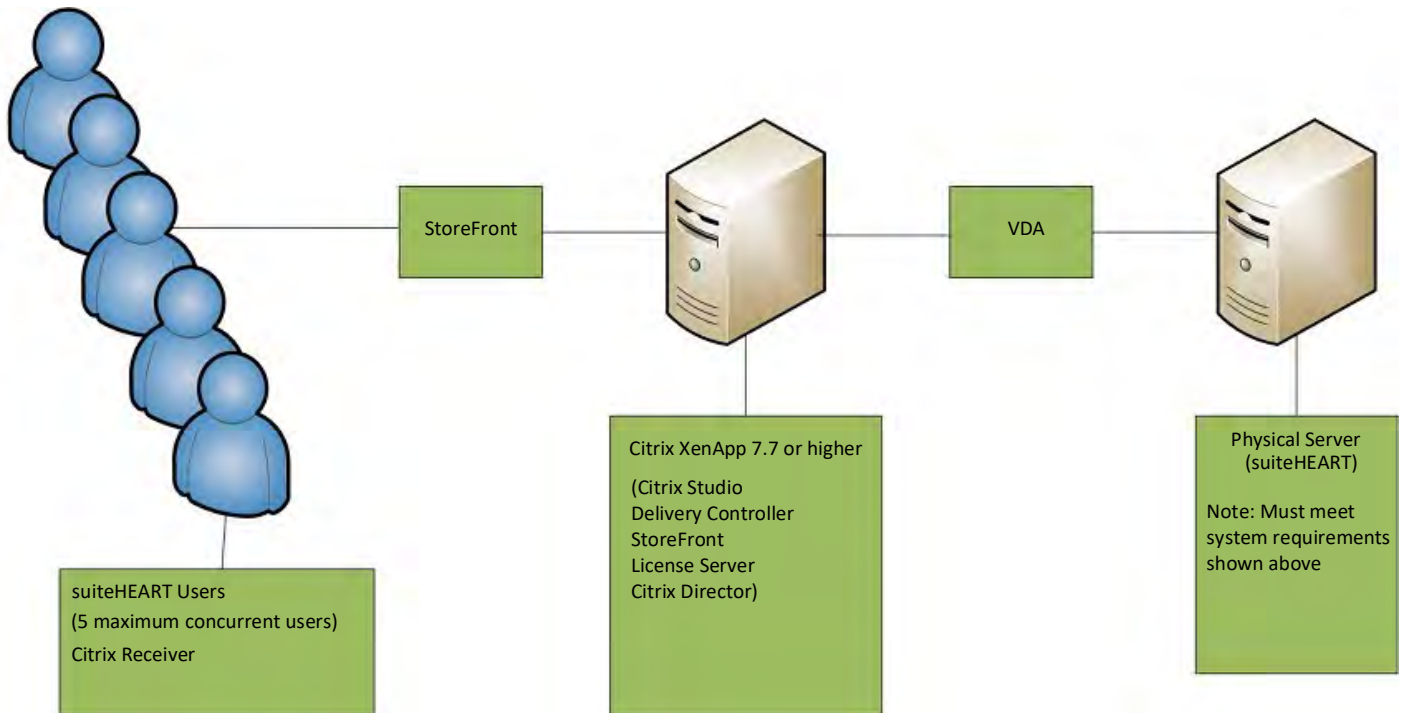Addendum

This addendum applies to suiteHEART® Software Instructions for Use and suiteDXT Instructions for Use.

suiteHEART® Software supports remote connectivity capability, as well as running in a multiuser environment. The remote connectivity capability allows suiteHEART® Software users remote access from computers not loaded with suiteHEART® Software. The remote user will have the ability to access the suiteHEART® Software system and have full functionality. Remote connectivity is supported for only a single concurrent remote connection, whereas the multiuser environment supports up to five concurrent users.

Remote Connectivity - Single User

suiteHEART® remote connectivity needs to be installed/configured/supported by your Information Technology personnel. Remote connectivity has been verified by NeoSoft using the following platforms:

- TeamViewer (version 14)
- Windows Remote Desktop (Microsoft Windows 7/10 Professional/Enterprise)

NOTE: When using Windows Remote Desktop, screen resolution on the remote computer must be set to 1920x1080 or higher video display resolution.

NOTE: When remotely accessing suiteHEART® Software using Windows Remote Desktop, the suiteHEART® 3D/4D viewer is not operational. The workaround is to use one of the other qualified remote access software package (TeamViewer) or an NVIDIA Quadro P4000/5000 graphics card.

The performance of the suiteHEART® Software via the remote connection cannot be guaranteed. Performance is dependent upon factors that are outside the scope of the suiteHEART® Software application. These factors include:

- Remote computer performance
- Internet download/upload speeds (Mbps)
- Network bandwidth speed (wired ethernet or wireless connection)

Multiuser Environment

The suiteDXT/suiteHEART® Software application allows up to five concurrent users when installed in a Citrix XenApp version 7.7 or higher environment.

Using suiteDXT Software in the Multiuser Environment

For complete installation instructions review the Installation Manual.

Configuration

Local Storage (i.e. Restart Service, Browse, Synchronize DB, Rebuild Database, and Permissions) configuration changes affect all system users and therefore should be managed by the system administrator for your facility.

Preprocessing


Preprocessing and Virtual Fellow™ configuration changes affect all system users and therefore should be managed by the System Administrator for your facility.

Launching a Study from suiteDXT

If a study has been opened in suiteHEART® Software by one user, that study is locked and a warning message will appear when a second user attempts to open the same study. The study can only be opened after it's closed by the first user. (Fig 1)

NOTE: The suiteHEART® Software will close after 60 minutes of inactivity and all analysis performed on open studies will be saved. The automatic closing of the suiteHEART® Software application will unlock open studies. It is recommended to always close the study after the analysis session is completed.

FIGURE 1. Launch Error

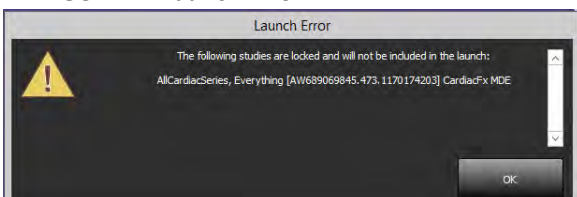
Query/Retrieve from suiteDXT

During Query/Retrieve the Network status indicator, in the lower right of the main DXT screen, will not be active (Fig. 2). The status of a Query/Retrieve can be viewed on the Query/Retrieve Window.

Any new study that has been networked or imported into suiteDXT may not appear automatically in the directory listing. Click on the **Refresh** button, on the main DXT screen, to refresh the study and series listings (Fig. 2).

FIGURE 2. Query/Retrieve

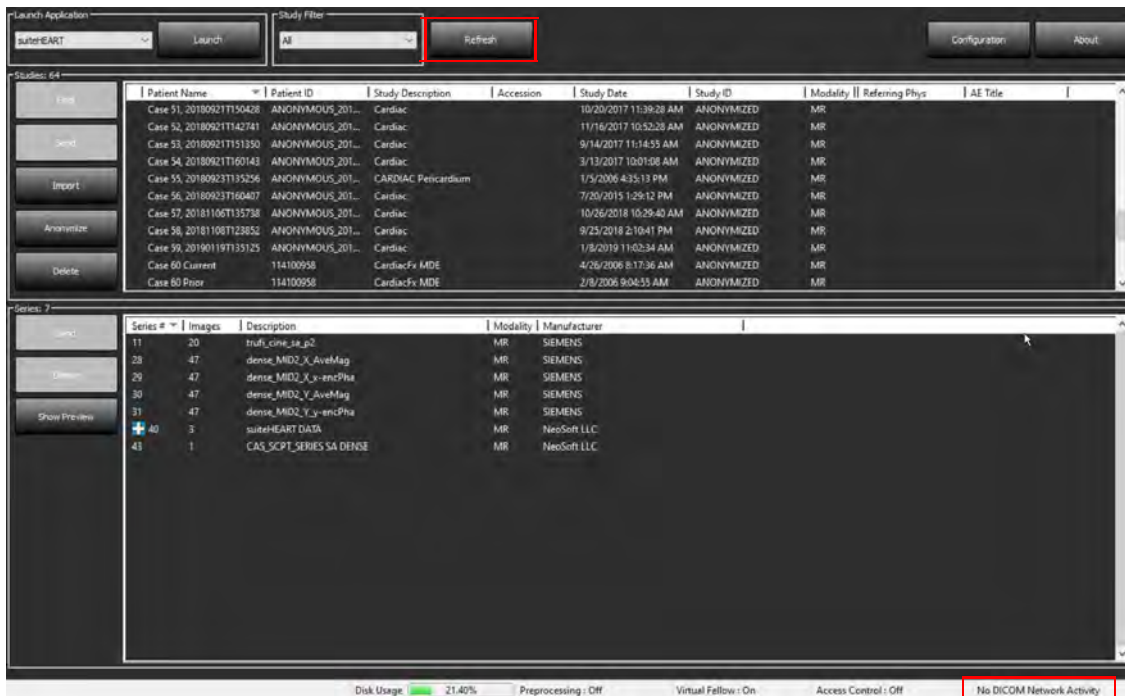
Deletion of a Study or Series from suiteDXT

Deletion of a study and its associated series that is opened and locked by another user is prohibited. A warning message (Fig 3) will appear if this is attempted. Deletion can only be performed when the study is closed by the other user.

NOTE: The suiteHEART® Software will close after 60 minutes of inactivity and all analysis performed on open studies will be saved. The automatic closing of the suiteHEART® Software application will unlock open studies. It is recommended to always close the study after the analysis session is completed.

FIGURE 3. Warning Message

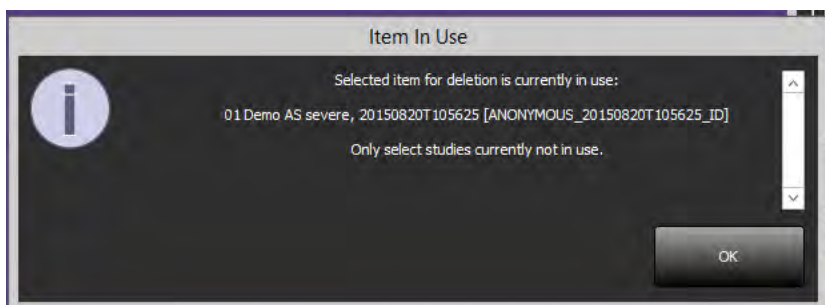Using suiteHEART® Software in the Multiuser Environment

For complete installation instructions review the Installation Manual.

Setting up Preferences

Editing report preferences and templates should be managed by the system administrator for your facility. Any changes to report preferences, templates, polar plot colors, print preferences and add measurements for reporting will impact all users and should only be performed by the system administrator.

Image Management Tools

When exporting file types such as Matlab, excel, JPEG, etcetera, verify the selection of the desired destination directory is on your local machine.

Report Database

Search queries created in the Report Database can be overwritten if other users are simultaneously creating search queries.

A fully qualified share folder (i.e. //192.168.3.82/test) must be used when defining a (SECONDARY_REPORT_FOLDER=) to be used to send approved exams from suiteHEART® Software in a multiuser environment.

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The screenshot shows the suiteHEART software interface. At the top, there is a 'Launch Application' section with a dropdown menu set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' section with a search box containing 'A' and a 'Refresh' button. Below this is a table of studies with columns for Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, Modality, Referring Phys, and AE Title. The table lists several cases, including Case 51 through Case 60. Below the studies table is a 'Series' section with a table of series with columns for Series #, Images, Description, Modality, and Manufacturer. The series table lists series 11 through 43. At the bottom of the interface, there is a status bar with indicators for Disk Usage (21.40%), Preprocessing (Off), Virtual Follow (On), Access Control (Off), and No DICOM Network Activity.

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The screenshot shows a dialog box titled 'Item In Use'. It contains an information icon on the left. The text inside the dialog box reads: 'Selected item for deletion is currently in use:' followed by '01 Demo AS severe, 20150820T105625 [ANONYMOUS_20150820T105625_ID]'. Below this, it says 'Only select studies currently not in use.' There is an 'OK' button at the bottom right of the dialog box.

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
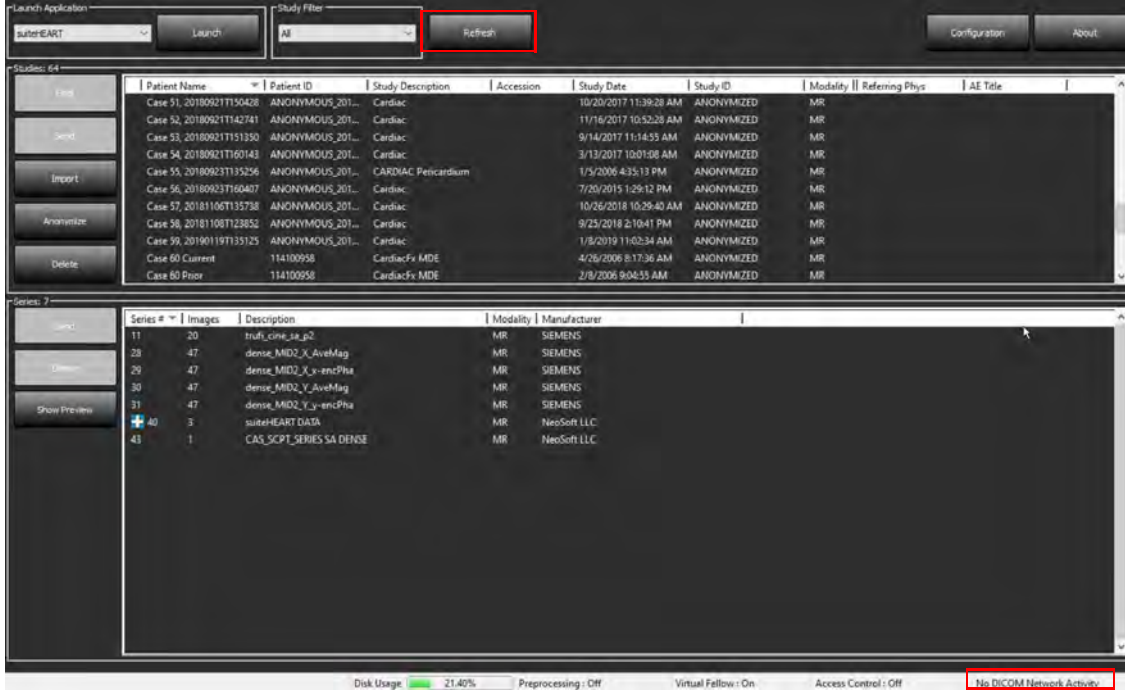
The screenshot shows the suiteDXT software interface. At the top, there is a 'Launch Application' section with a dropdown menu set to 'suiteHEART' and a 'Launch' button. To the right, there is a 'Study Filter' section with a dropdown menu set to 'All' and a 'Refresh' button highlighted with a red box. Below this, there is a table of studies with columns for Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, Modality, Referring Phys, and AE Title. The table contains 10 rows of data. Below the table, there is a 'Series 7' section with a table of series with columns for Series #, Images, Description, Modality, and Manufacturer. The table contains 5 rows of data. At the bottom of the interface, there is a status bar with 'Disk Usage: 21.40%', 'Preprocessing: Off', 'Virtual Follow: On', 'Access Control: Off', and 'No DICOM Network Activity' highlighted with a red box.

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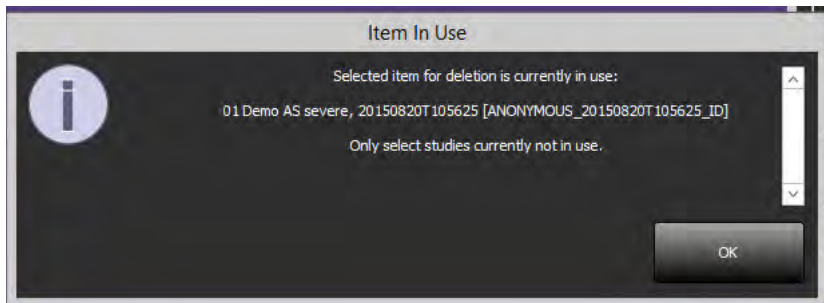The screenshot shows a dialog box titled 'Item In Use'. It contains an information icon on the left. The text inside the dialog box reads: 'Selected item for deletion is currently in use: 01 Demo AS severe, 20150820T105625 [ANONYMOUS_20150820T105625_ID]'. Below this, it says 'Only select studies currently not in use.' and there is an 'OK' button at the bottom right.

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