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NeoSoft, LLC

NEOSOFT

NS-03-039-0014 Rev 3
Copyright 2020 NeoSoft, LLC
All rights reserved

suiteHEART® Software

Addendum to Instructions for Use

NeoSoft, LLC

The logo for NeoSoft, featuring the word "NEO" in a bold, cyan-colored sans-serif font, followed by "SOFT" in a grey-colored sans-serif font. The letters are closely spaced and have a slight shadow effect.

NS-03-039-0014 Rev 3
Copyright 2020 NeoSoft, LLC
All rights reserved

Revision History

Rev	Date	Description of Change	Safety Related Update (Yes/No)
1	26APR2019	Updated for 5.0.0 product release. English manual previously NS-03-015-0006, Rev 7. Added French, German, Greek, and Italian languages which were previously: NS-03-015-0007 FR Rev 7 NS-03-015-0008 DE Rev 7 NS-03-038-0009 EL Rev 1 NS-03-038-0018 IT Rev 1	No
2	12NOV2019	Updated for 5.0.1 product release.	No
3	12OCT2020	Updated for 5.0.2 product release.	No

Manufacturer

NeoSoft, LLC
N27 W23910A Paul Road
Pewaukee, WI 53072 USA

Phone: 262-522-6120
website: www.neosoftllc.com

Sales: orders@neosoftmedical.com
Service: service@neosoftmedical.com

The Instructions for Use for this device are provided electronically in Portable Document Format, (.pdf). A pdf viewer is required to view the Instructions for Use. A paper copy of the Instructions for Use can be provided upon request, at no cost, within 7 calendar days, by emailing service@neosoftmedical.com.

Instructions for Use may be accessed in the following ways:


1. After launching the application, click "Help" or "About" from the main screen. Select the "Instructions for Use" option. The Instructions for Use will open in a pdf viewer.
2. If the original installation package received from NeoSoft is available, open the zip file and navigate to the "Documentation" folder, then the "Instructions for Use" folder and double click on the Instructions for Use.pdf in your language, denoted by EN - English, FR - French, DE - German, EL - Greek, IT - Italian, LT - Lithuanian, ES - Spanish, SV - Swedish, TR - Turkish, RO - Romanian, NL - Dutch, ZH_CN - Simplified Chinese, PT-PT - Portuguese, HU - Hungarian.
3. Navigate to the folder where the application is installed. Locate the "Documentation" folder, open the "Instructions for Use" folder and double click on the Instructions for Use.pdf in your language, denoted by EN - English, FR - French, DE - German, EL - Greek, IT - Italian, LT - Lithuanian, ES - Spanish, SV - Swedish, TR - Turkish, RO - Romanian, NL - Dutch, ZH_CN - Simplified Chinese, PT-PT - Portuguese, HU - Hungarian.
4. Electronic copies of Instructions for Use are also available at www.neosoftllc.com/neosoft/product_manuals/.

Medical device directive

This product conforms with the requirements of council directive 93/42/EEC concerning medical devices when it bears the following CE mark of Conformity:

European Representatives:

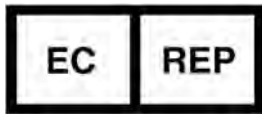EMERGO EUROPE
Prinsessegracht 20
2514 AP The Hague
The Netherlands

Canada:

Health Canada device license number: 99490

Malaysian Medical Device Registration Number: GB10979720-50888

Malaysia Authorized Representative:

Emergo Malaysia Sdn. Bhd.

Level 16, 1 Sentral Jalan Stesen Sentral 5 KL Sentral, 50470

Kuala Lumpur MALAYSIA


CAUTION: *United States federal law restricts this device to sale by, or on the order of, a physician.*

System Requirements

*Recommended video card(s); NVIDIA Quadro P4000 or P5000

suiteHEART® Software Minimum System Requirements (Multiuser Environment)
Supported Operating Systems: <ul style="list-style-type: none"> • Windows 2012 R2 Standard, Windows Server 2016 Standard, Edition (64 bit versions) or Windows Server 2019 Standard
Citrix XenApp 7.7 or higher
Microsoft® .NET Framework 4.0
Video card*(NVIDIA); 8GB RAM (minimum), minimum of 1664 CUDA Cores, support for OpenGL 4.0 and support for 1920 x 1080 resolution or higher
Monitor shall support 1920 x 1080 resolution or higher
Anti-Virus/Malware software installed
Minimum 64GB memory
Available hard drive space minimum 500GB
PDF viewer, Adobe Reader 11.0 or higher
An open network port for DICOM transfer over network
Static IP address (recommended)
Solid-State hard drive(s)
Intel® Core™ processor(s) (4 cores)
Citrix XenApp client displays must not exceed scaling/zoom levels greater than 125%

Citrix Installation (Example)

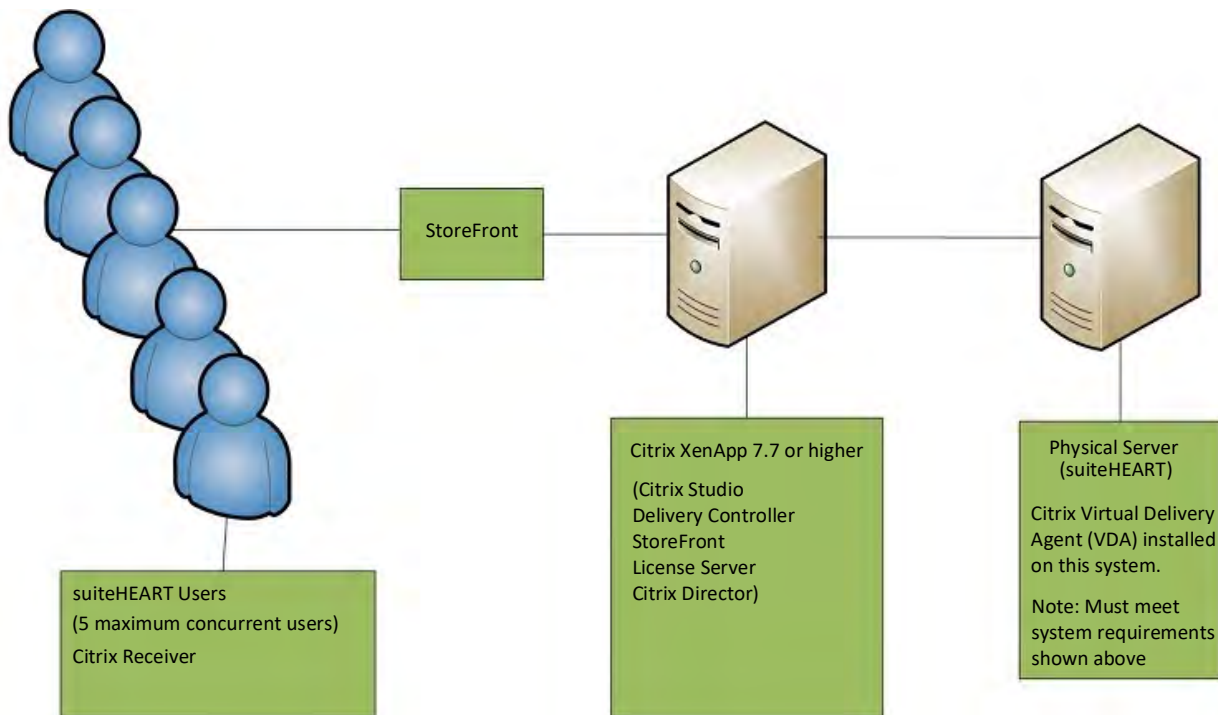Addendum

This addendum applies to suiteHEART® Software Instructions for Use and suiteDXT Instructions for Use.

suiteHEART® Software supports remote connectivity capability, as well as running in a multiuser environment. The remote connectivity capability allows suiteHEART® Software users remote access from computers not loaded with suiteHEART® Software. The remote user will have the ability to access the suiteHEART® Software system and have full functionality. Remote connectivity is supported for only a single concurrent remote connection, whereas the multiuser environment supports up to five concurrent users.

Remote Connectivity - Single User

suiteHEART® remote connectivity needs to be installed/configured/supported by your Information Technology personnel. Remote connectivity has been verified by NeoSoft using the following platforms:

- TeamViewer (version 15)
- Windows Remote Desktop (Microsoft Windows 7/10 Professional/Enterprise)

NOTE: When using Windows Remote Desktop, screen resolution on the remote computer must be set to 1920x1080 or higher video display resolution.

NOTE: When remotely accessing suiteHEART® Software using Windows Remote Desktop, the suiteHEART® 3D/4D viewer is not operational. The workaround is to use one of the other qualified remote access software package (TeamViewer) or an NVIDIA Quadro P4000/5000 graphics card.

The performance of the suiteHEART® Software via the remote connection cannot be guaranteed. Performance is dependent upon factors that are outside the scope of the suiteHEART® Software application. These factors include:

- Remote computer performance
- Internet download/upload speeds (Mbps)
- Network bandwidth speed (wired ethernet or wireless connection)

Multiuser Environment

The suiteDXT/suiteHEART® Software application allows up to five concurrent users when installed in a Citrix XenApp version 7.7 or higher environment.

Using suiteDXT Software in the Multiuser Environment

For complete installation instructions review the Installation Manual.

Configuration

Local Storage (i.e. Restart Service, Browse, Synchronize DB and Permissions) configuration changes affect all system users and therefore should be managed by the system administrator for your facility.

Preprocessing


Preprocessing and Virtual Fellow™ configuration changes affect all system users and therefore should be managed by the System Administrator for your facility.

Launching a Study from suiteDXT

If a study has been opened in suiteHEART® Software by one user, that study is locked and a warning message will appear when a second user attempts to open the same study. The study can only be opened after it's closed by the first user unless Read Only mode is enabled. In this mode, subsequent launches will be opened but in Read Only format. (Fig 1)

NOTE: The suiteHEART® Software will close after 60 minutes of inactivity and all analysis performed on open studies will be saved. The automatic closing of the suiteHEART® Software application will unlock open studies. It is recommended to always close the study after the analysis session is completed.

FIGURE 1. Launch Error

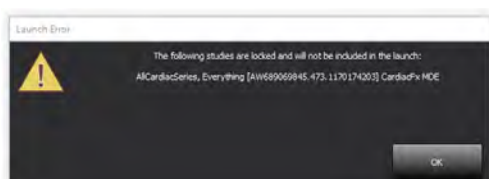Query/Retrieve from suiteDXT

During Query/Retrieve the Network status indicator, in the lower right of the main DXT screen, will not be active (Fig. 2). The status of a Query/Retrieve can be viewed on the Query/Retrieve Window.

Any new study that has been networked or imported into suiteDXT may not appear automatically in the directory listing. Click on the **Refresh** button, on the main DXT screen, to refresh the study and series listings (Fig. 2).

FIGURE 2. Query/Retrieve

The screenshot shows the suiteDXT interface. At the top, there is a 'Launch Application' dropdown set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' dropdown set to 'All' and a 'Refresh' button, which is highlighted with a red box. Further right are 'Configuration' and 'About' buttons. Below this is a section for 'Studies: 52' with a search bar and a table of study details. Below that is a section for 'Series: 24' with a search bar and a table of series details. At the bottom, there is a status bar showing 'Disk Usage 67.15%', 'Preprocessing: On', 'Virtual Fellow®: On', and 'No DICOM Network Activity'.

Patient Name	Patient ID	Study Description	Accession	Study Date	Study ID	Modality
Case 37 Current	7878	Cardiac		2/22/2017 3:04:56 PM	ANONYMIZED	MR
Case 37 Prior	7878	Cardiac		2/14/2013 12:00:00 AM	ANONYMIZED	MR
Case 43, 20180922T110702	ANONYMOUS_20180922T110702_ID	Cardiac		9/6/2017 12:23:09 PM	ANONYMIZED	MR
Case 44, 20180922T153755	ANONYMOUS_20180922T153755_ID	Cardiac		11/29/2017 9:36:20 AM	ANONYMIZED	MR
Case 45, 20180922T180222	ANONYMOUS_20180922T180222_ID	Cardiac		10/17/2017 11:27:47 AM	ANONYMIZED	MR
Case 46, 20180922T171615	ANONYMOUS_20180922T171615_ID	Cardiac/MRA: Chest		9/7/2018 10:26:26 AM	ANONYMIZED	MR
Case 47, 20180922T172103	ANONYMOUS_20180922T172103_ID	Cardiac		7/28/2017 5:00:02 PM	ANONYMIZED	MR
Case 48	ANONYMOUS_20191008T163809_ID	Cardiac		4/18/2019 9:56:52 AM	ANONYMIZED	MR
Case 49	ANONYMOUS_20191009T162750_ID	Cardiac		8/1/2016 9:42:03 AM	ANONYMIZED	MR
Case 50	ANONYMOUS_20191010T111927_ID	Cardiac UVA BILCHL...		1/1/2022 12:00:00 AM	ANONYMIZED	MR
Case 51, 20180921T150428	ANONYMOUS_20180921T150428_ID	Cardiac		10/20/2017 11:39:28 AM	ANONYMIZED	MR


Series #	Images	Description	Modality	Manufacturer	Transferred To	Transferred Date
3003	2400	4DFLOW - CE - Wholechest - SI Flow	MR	GE MEDICAL SYSTEMS		
3002	2400	4DFLOW - CE - Wholechest - AP Flow	MR	GE MEDICAL SYSTEMS		
3001	2400	4DFLOW - CE - Wholechest - LR Flow	MR	GE MEDICAL SYSTEMS		
3000	2400	4DFLOW - CE - Wholechest - Anatomy	MR	GE MEDICAL SYSTEMS		
38	60	15 Ao PC	MR	GE MEDICAL SYSTEMS		
37	60	14 PA PC	MR	GE MEDICAL SYSTEMS		
36	60	13 Ao(BCT) PC	MR	GE MEDICAL SYSTEMS		
34	1	PS:2D PSMDE	MR	GE MEDICAL SYSTEMS		
30	480	4DFLOW - CE - Wholechest - Preview	MR	GE MEDICAL SYSTEMS		

Deletion of a Study or Series from suiteDXT

Deletion of a study and its associated series that is opened and locked by another user is prohibited. A warning message (Fig 3) will appear if this is attempted. Deletion can only be performed when the study is closed by the other user.

NOTE: The suiteHEART® Software will close after 60 minutes of inactivity and all analysis performed on open studies will be saved. The automatic closing of the suiteHEART® Software application will unlock open studies. It is recommended to always close the study after the analysis session is completed.

FIGURE 3. Warning Message

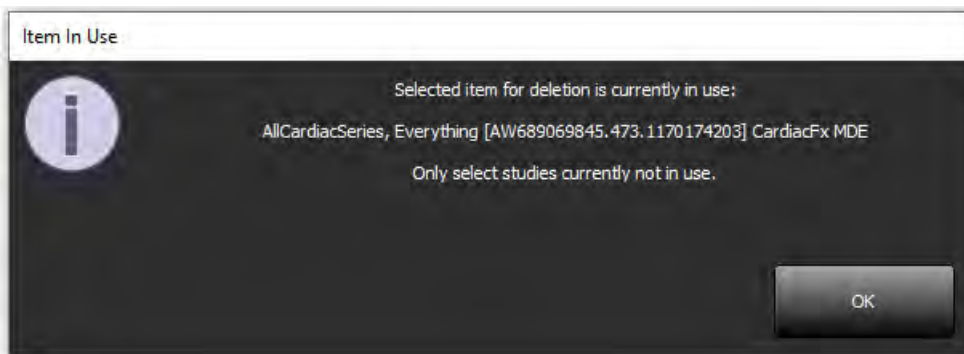Using suiteHEART® Software in the Multiuser Environment

For complete installation instructions review the Installation Manual.

Setting up Preferences

Editing preferences should be managed by the system administrator for your facility.

Image Management Tools

When exporting file types such as Matlab, excel, JPEG, etcetera, verify the selection of the desired destination directory is on your local machine.

Report Database

Search queries created in the Report Database can be overwritten if other users are simultaneously creating search queries.

A fully qualified share folder (i.e. //192.168.3.82/test) must be used when defining a (SECONDARY_REPORT_FOLDER=) to be used to send approved exams from suiteHEART® Software in a multiuser environment.

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
FIGURE 2. Interrogation/Récupération

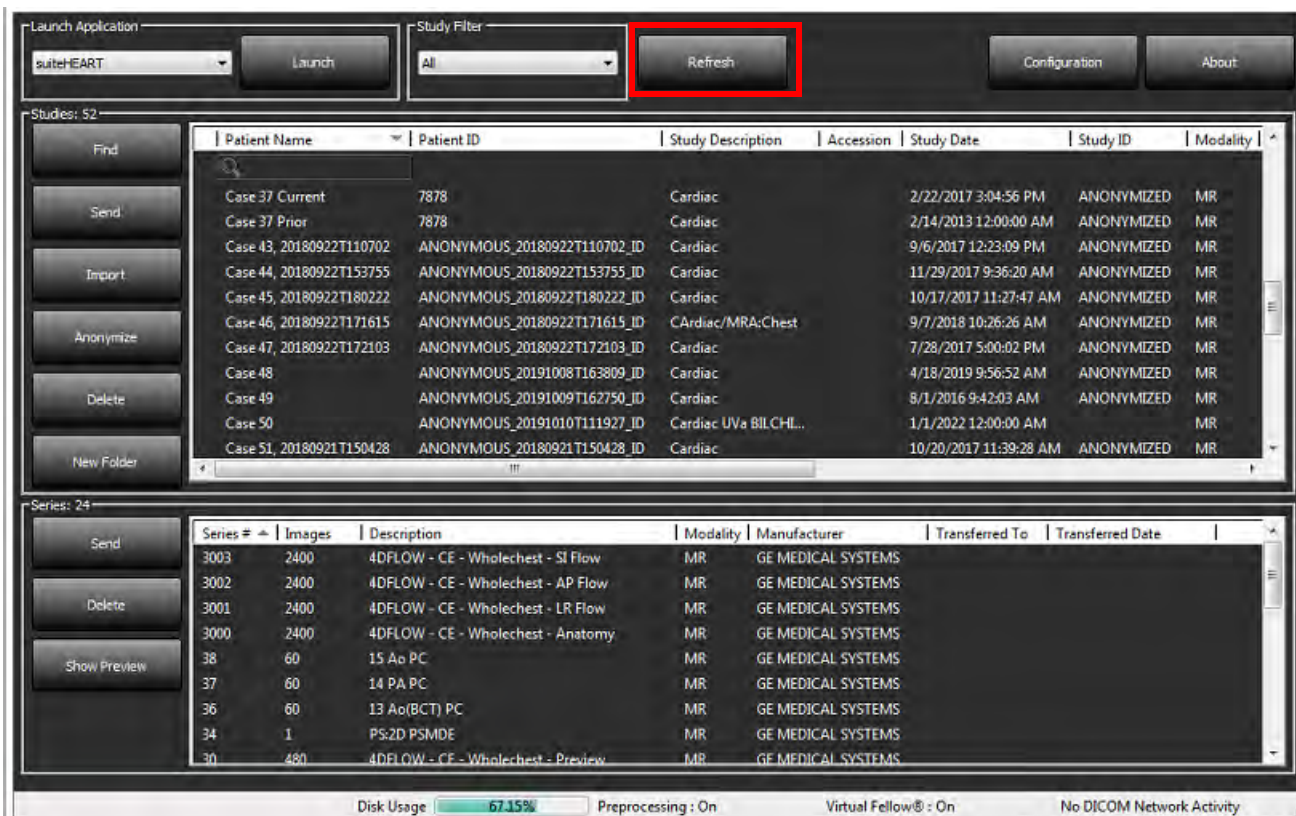The screenshot shows the suiteHEART application interface. At the top, there is a 'Launch Application' section with a dropdown menu set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' dropdown set to 'All' and a 'Refresh' button, which is highlighted with a red rectangle. Further right are 'Configuration' and 'About' buttons. Below this is a 'Studies: 52' section with a search bar and a list of studies. The list has columns for Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, and Modality. Below the studies list is a 'Series: 24' section with a search bar and a list of series. The series list has columns for Series #, Images, Description, Modality, Manufacturer, Transferred To, and Transferred Date. At the bottom of the interface, there is a status bar showing 'Disk Usage 67.15%', 'Preprocessing : On', 'Virtual Fellow® : On', and 'No DICOM Network Activity'.

Patient Name	Patient ID	Study Description	Accession	Study Date	Study ID	Modality
Case 37 Current	7878	Cardiac		2/22/2017 3:04:56 PM	ANONYMIZED	MR
Case 37 Prior	7878	Cardiac		2/14/2013 12:00:00 AM	ANONYMIZED	MR
Case 43, 20180922T110702	ANONYMOUS_20180922T110702_ID	Cardiac		9/6/2017 12:23:09 PM	ANONYMIZED	MR
Case 44, 20180922T153755	ANONYMOUS_20180922T153755_ID	Cardiac		11/29/2017 9:36:20 AM	ANONYMIZED	MR
Case 45, 20180922T180222	ANONYMOUS_20180922T180222_ID	Cardiac		10/17/2017 11:27:47 AM	ANONYMIZED	MR
Case 46, 20180922T171615	ANONYMOUS_20180922T171615_ID	CARDIAC/MRA;Chest		9/7/2018 10:26:26 AM	ANONYMIZED	MR
Case 47, 20180922T172103	ANONYMOUS_20180922T172103_ID	Cardiac		7/28/2017 5:00:02 PM	ANONYMIZED	MR
Case 48	ANONYMOUS_20191008T163809_ID	Cardiac		4/18/2019 9:56:52 AM	ANONYMIZED	MR
Case 49	ANONYMOUS_20191009T162750_ID	Cardiac		8/1/2016 9:42:03 AM	ANONYMIZED	MR
Case 50	ANONYMOUS_20191010T111927_ID	Cardiac Uva BILCHL...		1/1/2022 12:00:00 AM		MR
Case 51, 20180921T150428	ANONYMOUS_20180921T150428_ID	Cardiac		10/20/2017 11:39:28 AM	ANONYMIZED	MR

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
FIGURA 2. Ricerca/recupero

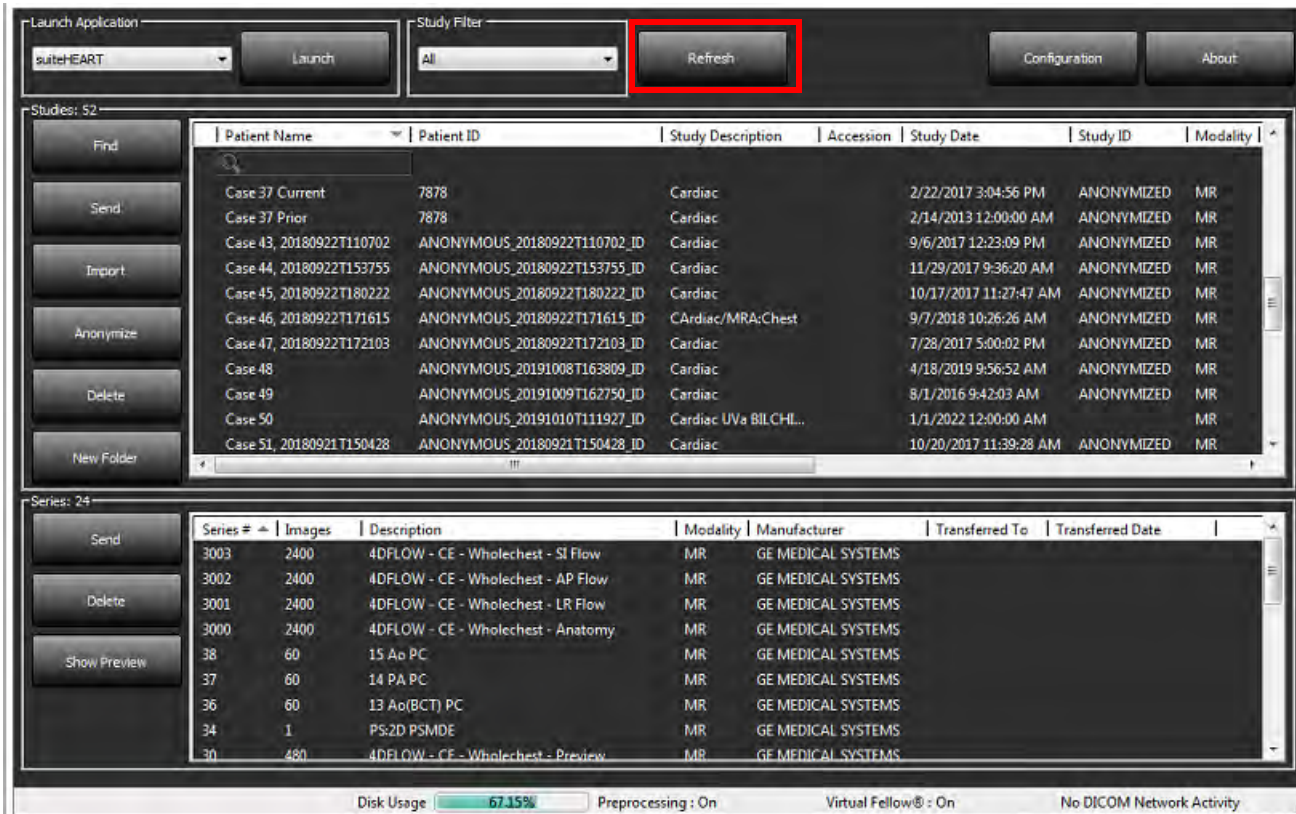The screenshot shows the suiteDXT interface with the 'Refresh' button highlighted in red. The interface displays a list of studies and series.

Studies: 52

Patient Name	Patient ID	Study Description	Accession	Study Date	Study ID	Modality
Case 37 Current	7878	Cardiac		2/22/2017 3:04:56 PM	ANONYMIZED	MR
Case 37 Prior	7878	Cardiac		2/14/2013 12:00:00 AM	ANONYMIZED	MR
Case 43, 20180922T110702	ANONYMOUS_20180922T110702_ID	Cardiac		9/6/2017 12:23:09 PM	ANONYMIZED	MR
Case 44, 20180922T153755	ANONYMOUS_20180922T153755_ID	Cardiac		11/29/2017 9:36:20 AM	ANONYMIZED	MR
Case 45, 20180922T180222	ANONYMOUS_20180922T180222_ID	Cardiac		10/17/2017 11:27:47 AM	ANONYMIZED	MR
Case 46, 20180922T171615	ANONYMOUS_20180922T171615_ID	Cardiac/MRA:Chest		9/7/2018 10:26:26 AM	ANONYMIZED	MR
Case 47, 20180922T172103	ANONYMOUS_20180922T172103_ID	Cardiac		7/28/2017 5:00:02 PM	ANONYMIZED	MR
Case 48	ANONYMOUS_20191008T163809_ID	Cardiac		4/18/2019 9:56:52 AM	ANONYMIZED	MR
Case 49	ANONYMOUS_20191009T162750_ID	Cardiac		8/1/2016 9:42:03 AM	ANONYMIZED	MR
Case 50	ANONYMOUS_20191010T111927_ID	Cardiac UVa BILCHL...		1/1/2022 12:00:00 AM		MR
Case 51, 20180921T150428	ANONYMOUS_20180921T150428_ID	Cardiac		10/20/2017 11:39:28 AM	ANONYMIZED	MR

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Disk Usage: 67.15% | Preprocessing: On | Virtual Fellow®: On | No DICOM Network Activity

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The screenshot shows the suiteDXT interface. At the top, there is a 'Launch Application' section with a dropdown menu set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' section with a dropdown menu set to 'All' and a 'Refresh' button highlighted with a red box. Further right are 'Configuration' and 'About' buttons. Below this is a 'Studies: 52' section with a search bar and a table of study records. The table has columns for Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, and Modality. Below the studies section is a 'Series: 24' section with a search bar and a table of series records. The table has columns for Series #, Images, Description, Modality, Manufacturer, Transferred To, and Transferred Date. At the bottom of the interface, there is a status bar showing 'Disk Usage: 67.15%', 'Preprocessing: On', 'Virtual Fellow®: On', and 'No DICOM Network Activity'.

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
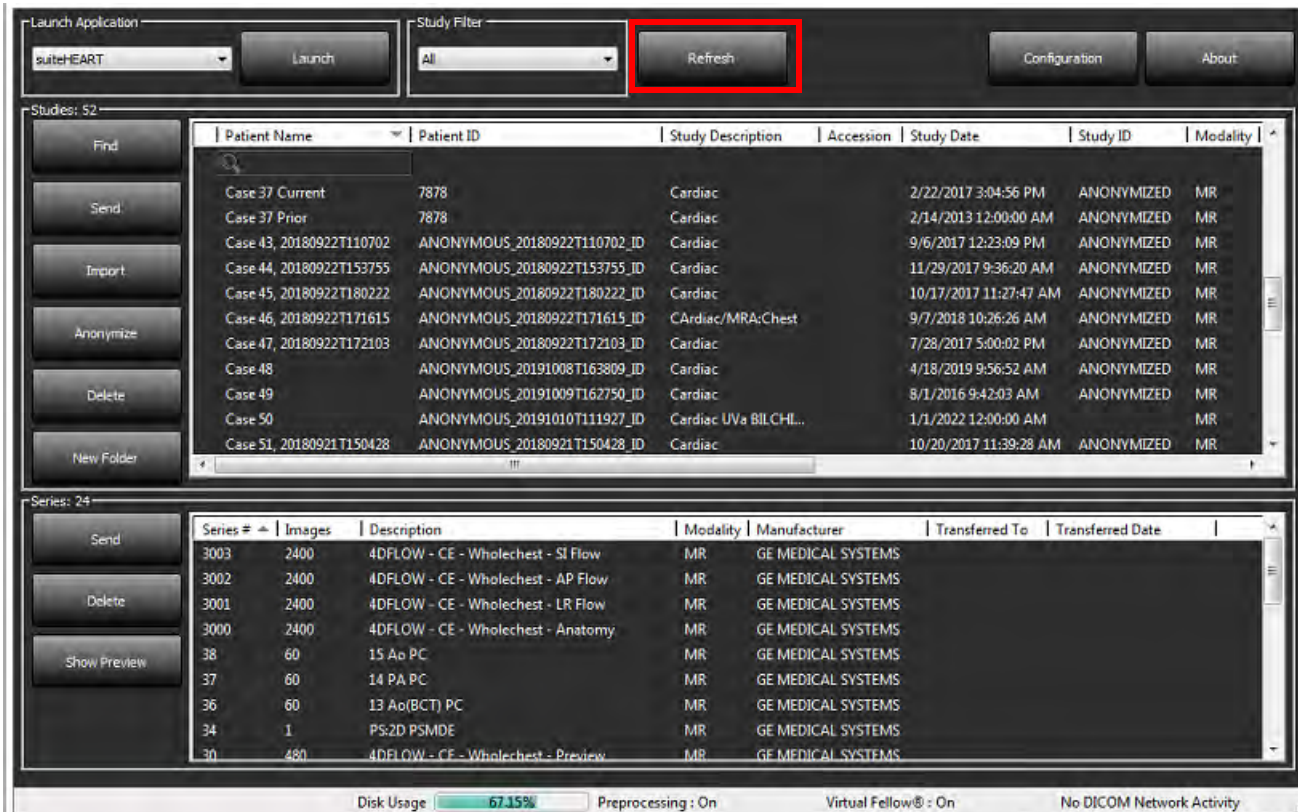
The screenshot shows the suiteDXT interface. At the top, there is a 'Launch Application' dropdown set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' dropdown set to 'All' and a 'Refresh' button highlighted with a red box. Further right are 'Configuration' and 'About' buttons. Below this is a table of studies with columns: Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, and Modality. The table lists 11 cases, including 'Case 37 Current', 'Case 37 Prior', and various 'Case 43' through 'Case 51' entries. Below the studies table is a table of series with columns: Series #, Images, Description, Modality, Manufacturer, Transferred To, and Transferred Date. It lists 10 series, including '3003 2400 4DFLOW - CE - Wholechest - SI Flow' and '3002 2400 4DFLOW - CE - Wholechest - AP Flow'. At the bottom, there is a status bar showing 'Disk Usage 67.15%', 'Preprocessing : On', 'Virtual Fellow® : On', and 'No DICOM Network Activity'.

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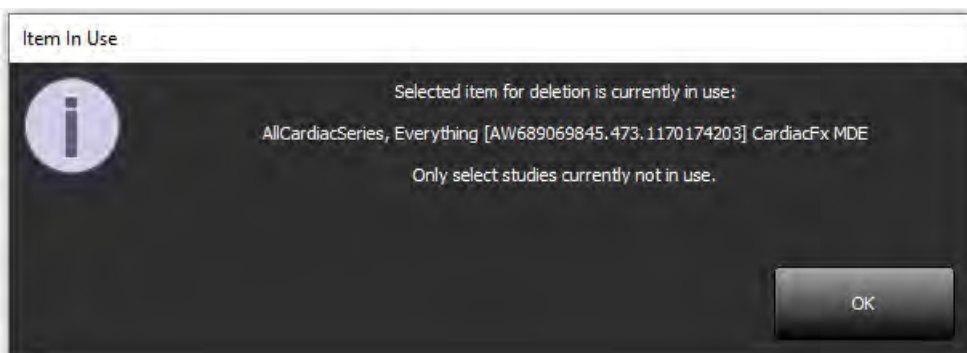The screenshot shows an error dialog box titled 'Item In Use'. It contains an information icon (i) and the following text: 'Selected item for deletion is currently in use: AllCardiacSeries, Everything [AW689069845.473.1170174203] CardiacFx MDE'. Below this, it says 'Only select studies currently not in use.' and there is an 'OK' button at the bottom right.

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
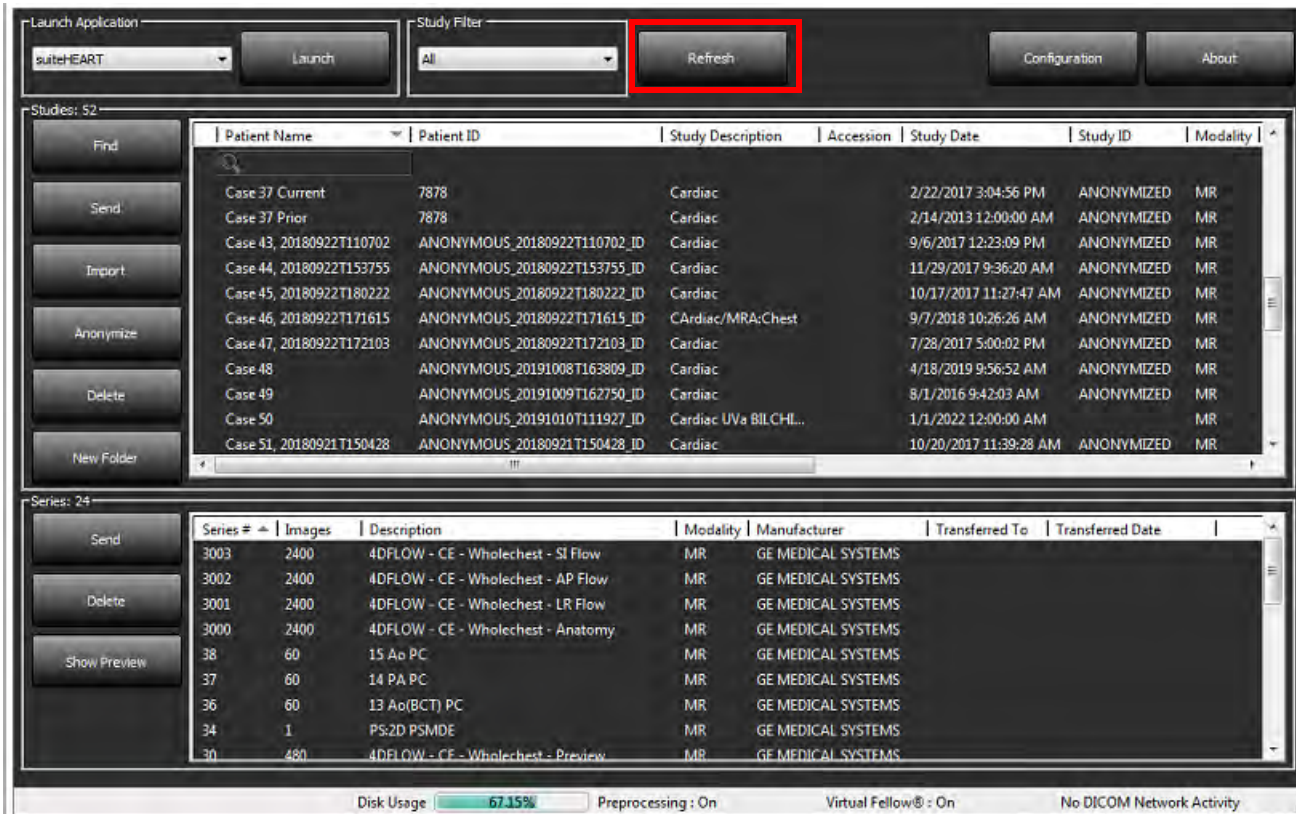
The screenshot shows the suiteDXT interface. At the top, there is a 'Launch Application' dropdown set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' dropdown set to 'All' and a 'Refresh' button, which is highlighted with a red box. Further right are 'Configuration' and 'About' buttons. Below this is a 'Studies: 52' section with a search bar and a list of studies. The list has columns for Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, and Modality. Below the studies list is a 'Series: 24' section with a search bar and a list of series. The list has columns for Series #, Images, Description, Modality, Manufacturer, Transferred To, and Transferred Date. At the bottom of the interface, there is a status bar showing 'Disk Usage 67.15%', 'Preprocessing : On', 'Virtual Fellow® : On', and 'No DICOM Network Activity'.

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BILD 3. Varningsmeddelande

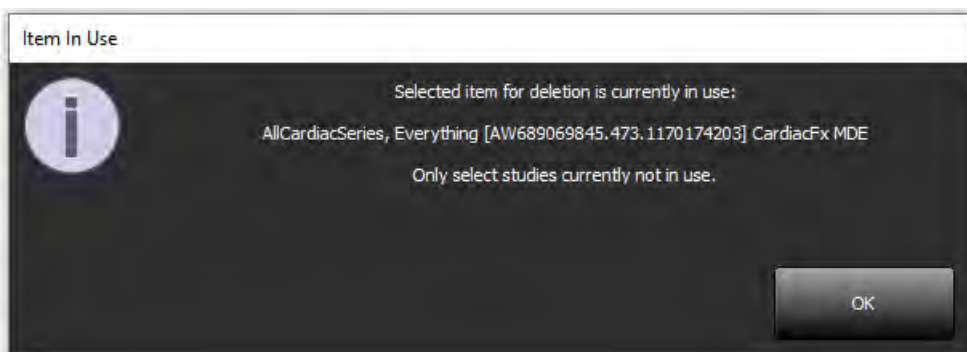The screenshot shows a warning dialog box titled 'Item In Use'. It contains an information icon on the left. The text inside the dialog reads: 'Selected item for deletion is currently in use: AllCardiacSeries, Everything [AW689069845.473.1170174203] CardiacFx MDE. Only select studies currently not in use.' There is an 'OK' button at the bottom right of the dialog.

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
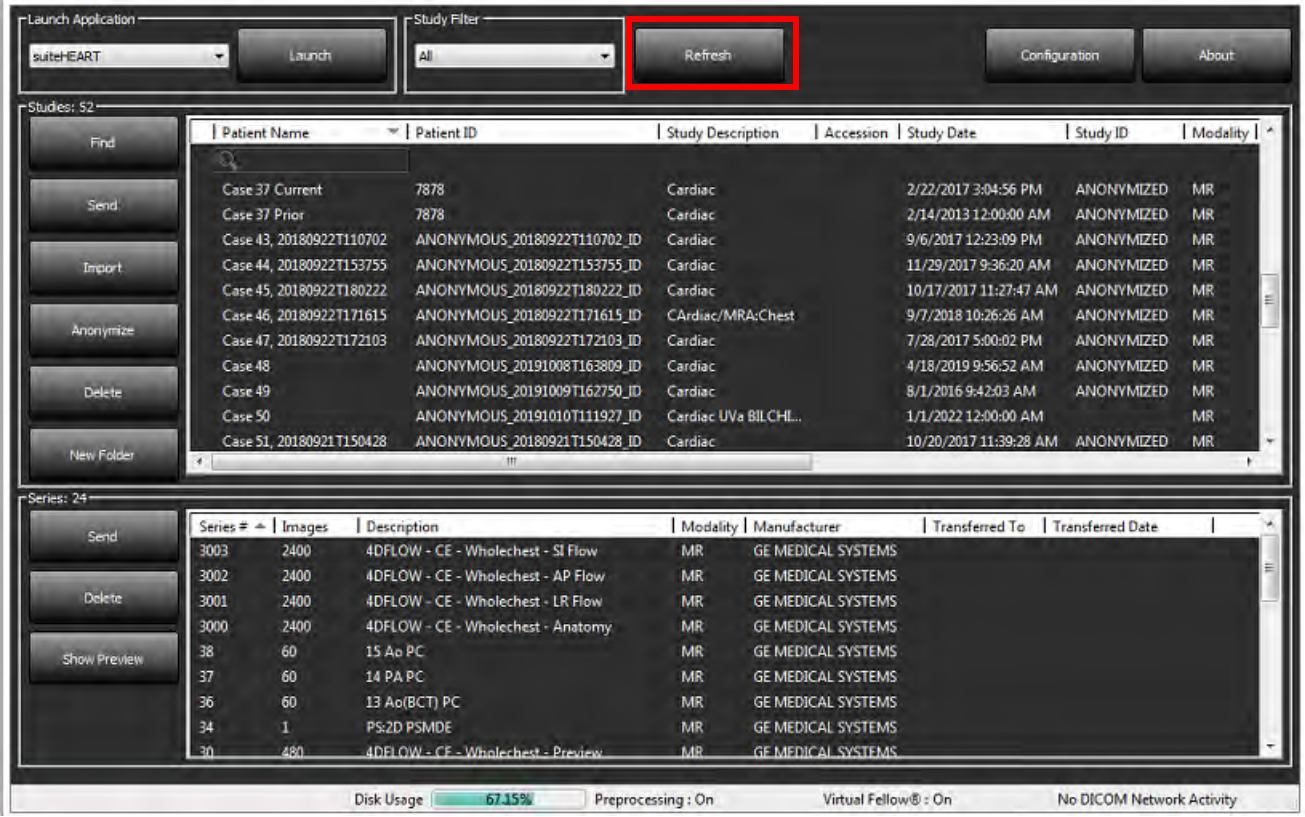
The screenshot shows the suiteDXT software interface. At the top, there is a 'Launch Application' section with a dropdown menu set to 'suiteHEART' and a 'Launch' button. To the right, there is a 'Study Filter' dropdown set to 'All' and a 'Refresh' button, which is highlighted with a red box. Further right are 'Configuration' and 'About' buttons. Below this, there is a 'Studies: 52' section with a search bar and a list of studies. The list has columns for Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, and Modality. Below the studies list is a 'Series: 24' section with a search bar and a list of series. The list has columns for Series #, Images, Description, Modality, Manufacturer, Transferred To, and Transferred Date. At the bottom of the interface, there is a status bar showing 'Disk Usage 67.15%', 'Preprocessing : On', 'Virtual Fellow® : On', and 'No DICOM Network Activity'.

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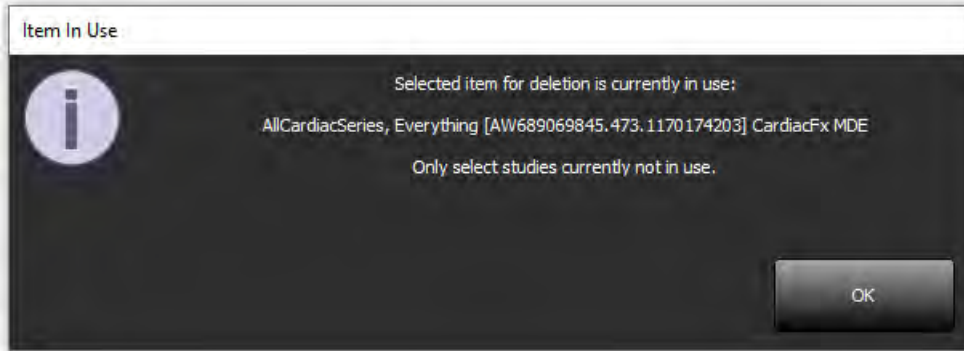The screenshot shows a warning message dialog box titled 'Item In Use'. The message text is: 'Selected item for deletion is currently in use: AllCardiacSeries, Everything [AW689069845.473, 1170174203] CardiacFx MDE. Only select studies currently not in use.' There is an 'OK' button at the bottom right of the dialog box.

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FIGUUR 2. Opvragen/Ophalen

The screenshot shows the suiteHEART software interface. At the top, there is a 'Launch Application' dropdown set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' dropdown set to 'All' and a 'Refresh' button, which is highlighted with a red box. Further right are 'Configuration' and 'About' buttons. Below this is a 'Studies: 52' section with a table of study data. Below that is a 'Series: 24' section with a table of series data. At the bottom, there is a status bar showing 'Disk Usage 67.15%', 'Preprocessing : On', 'Virtual Fellow® : On', and 'No DICOM Network Activity'.

Patient Name	Patient ID	Study Description	Accession	Study Date	Study ID	Modality
Case 37 Current	7878	Cardiac		2/22/2017 3:04:56 PM	ANONYMIZED	MR
Case 37 Prior	7878	Cardiac		2/14/2013 12:00:00 AM	ANONYMIZED	MR
Case 43, 20180922T110702	ANONYMOUS_20180922T110702_ID	Cardiac		9/6/2017 12:23:09 PM	ANONYMIZED	MR
Case 44, 20180922T153755	ANONYMOUS_20180922T153755_ID	Cardiac		11/29/2017 9:36:20 AM	ANONYMIZED	MR
Case 45, 20180922T180222	ANONYMOUS_20180922T180222_ID	Cardiac		10/17/2017 11:27:47 AM	ANONYMIZED	MR
Case 46, 20180922T171615	ANONYMOUS_20180922T171615_ID	CARDIAC/MRA:CHEST		9/7/2018 10:26:26 AM	ANONYMIZED	MR
Case 47, 20180922T172103	ANONYMOUS_20180922T172103_ID	Cardiac		7/28/2017 5:00:02 PM	ANONYMIZED	MR
Case 48	ANONYMOUS_20191008T163809_ID	Cardiac		4/18/2019 9:56:52 AM	ANONYMIZED	MR
Case 49	ANONYMOUS_20191009T162750_ID	Cardiac		8/1/2016 9:42:03 AM	ANONYMIZED	MR
Case 50	ANONYMOUS_20191010T111927_ID	Cardiac UVA BILCHL...		1/1/2022 12:00:00 AM	ANONYMIZED	MR
Case 51, 20180921T150428	ANONYMOUS_20180921T150428_ID	Cardiac		10/20/2017 11:39:28 AM	ANONYMIZED	MR

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FIGURA 2. Consultar/Recuperar

Patient Name	Patient ID	Study Description	Accession	Study Date	Study ID	Modality
Case 37 Current	7878	Cardiac		2/22/2017 3:04:56 PM	ANONYMIZED	MR
Case 37 Prior	7878	Cardiac		2/14/2013 12:00:00 AM	ANONYMIZED	MR
Case 43, 20180922T110702	ANONYMOUS_20180922T110702_ID	Cardiac		9/6/2017 12:23:09 PM	ANONYMIZED	MR
Case 44, 20180922T153755	ANONYMOUS_20180922T153755_ID	Cardiac		11/29/2017 9:36:20 AM	ANONYMIZED	MR
Case 45, 20180922T180222	ANONYMOUS_20180922T180222_ID	Cardiac		10/17/2017 11:27:47 AM	ANONYMIZED	MR
Case 46, 20180922T171615	ANONYMOUS_20180922T171615_ID	CArdiac/MRA:Chest		9/7/2018 10:26:26 AM	ANONYMIZED	MR
Case 47, 20180922T172103	ANONYMOUS_20180922T172103_ID	Cardiac		7/28/2017 5:00:02 PM	ANONYMIZED	MR
Case 48	ANONYMOUS_20191008T163809_ID	Cardiac		4/18/2019 9:56:52 AM	ANONYMIZED	MR
Case 49	ANONYMOUS_20191009T162750_ID	Cardiac		8/1/2016 9:42:03 AM	ANONYMIZED	MR
Case 50	ANONYMOUS_20191010T111927_ID	Cardiac UVA BILCHL...		1/1/2022 12:00:00 AM		MR
Case 51, 20180921T150428	ANONYMOUS_20180921T150428_ID	Cardiac		10/20/2017 11:39:28 AM	ANONYMIZED	MR

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The screenshot shows the suiteHEART software interface. At the top, there is a 'Launch Application' dropdown set to 'suiteHEART' and a 'Launch' button. To the right is a 'Study Filter' dropdown set to 'All' and a 'Refresh' button, which is highlighted with a red box. Further right are 'Configuration' and 'About' buttons. Below this is a 'Studies: 52' section with a search bar and buttons for 'Find', 'Send', 'Import', 'Anonymize', 'Delete', and 'New Folder'. The main area contains a table of studies with columns: Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, and Modality. Below that is a 'Series: 24' section with buttons for 'Send', 'Delete', and 'Show Preview'. It contains a table of series with columns: Series #, Images, Description, Modality, Manufacturer, Transferred To, and Transferred Date. At the bottom, there is a status bar showing 'Disk Usage: 67.15%', 'Preprocessing: On', 'Virtual Fellow: On', and 'No DICOM Network Activity'.

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