[suiteHEART® Software Installation Manual](#)

[suiteHEART® Software Installationshandbuch](#)

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Installation Manual

NeoSoft, LLC

NS-03-039-0013 Rev 4
Copyright 2021 NeoSoft, LLC
All rights reserved

Revision History

Rev	Date	Description of Change	Safety Related Update (Yes/No)
1	23APR2019	Updated for 5.0.0 product release. English manual previously NS-03-028-0003, Rev 4. Added French, German, Greek, and Italian languages which were previously: NS-03-028-0005 FR Rev 4 NS-03-028-0004 DE Rev 4 NS-03-038-0011 EL Rev 1 NS-03-038-0020 IT Rev 1	No
2	11NOV2019	Updated for 5.0.1 product release.	No
3	20OCT2020	Updated for 5.0.2 product release.	No
4	14JUNE2021	Added MD symbol, EU importer reference, incident reporting information.	No

Manufacturer

NeoSoft, LLC
N27 W23910A Paul Road
Pewaukee, WI 53072 USA

Phone: 262-522-6120
website: www.neosoftllc.com

Sales: orders@neosoftmedical.com
Service: service@neosoftmedical.com

The Instructions for Use for this device are provided electronically in Portable Document Format, (.pdf). A pdf viewer is required to view the Instructions for Use. A paper copy of the Instructions for Use can be provided upon request, at no cost, within 7 calendar days, by emailing service@neosoftmedical.com.

Instructions for Use may be accessed in the following ways:

1. After launching the application, click “Help” or “About” from the main screen. Select the “Instructions for Use” option. The Instructions for Use will open in a pdf viewer.
2. If the original installation package received from NeoSoft is available, open the zip file and navigate to the “Documentation” folder, then the “Instructions for Use” folder and double click on the Instructions for Use.pdf in your language, denoted by EN - English, FR - French, DE - German, EL - Greek, IT - Italian, LT - Lithuanian, ES - Spanish, SV - Swedish, TR - Turkish, RO - Romanian, NL - Dutch, PT-PT - Portuguese, HU - Hungarian.
3. Navigate to the folder where the application is installed. Locate the “Documentation” folder, open the “Instructions for Use” folder and double click on the Instructions for Use.pdf in your language, denoted by EN - English, FR - French, DE - German, EL - Greek, IT - Italian, LT - Lithuanian, ES - Spanish, SV - Swedish, TR - Turkish, RO - Romanian, NL - Dutch, PT-PT - Portuguese, HU - Hungarian.
4. Electronic copies of Instructions for Use are also available at www.neosoftllc.com/neosoft/product_manuals/.

Table of Contents

Safety	1
<i>Introduction</i>	1
<i>Indications for Use</i>	1
<i>Intended Use</i>	2
<i>Supported DICOM Image Formats</i>	2
<i>Safety Notices</i>	3
<i>Equipment Hazards</i>	3
suiteDXT Installation & Configuration	4
Prerequisites	4
Installation	5
Configuration	5
Local Storage Tab	8
<i>Auto Delete</i>	9
Application Launch Tab	9
Feature Licensing Tab	10
<i>Request a License</i>	11
<i>Apply the License</i>	11
Tablet Mode	11
suiteHEART® Software Installation	12
Prerequisites	12
Installation	13
<i>Software</i>	13
<i>Documentation</i>	13
Anti-Virus Exclusions	13
Application Inactivity Setting	14
Licensing Schemes	14
<i>Remote Connectivity</i>	15
Install Verification	15

Medical device directive

This product conforms with the requirements of council directive 93/42/EEC concerning medical devices when it bears the following CE mark of Conformity:

This product is a Medical Device:

European Representatives:

EMERGO EUROPE
Prinsessegracht 20
2514 AP The Hague
The Netherlands

EU Importer:

MedEnvoy Global B.V.
Pr. Margrietplantsoen 33, Suite 123,
2595 AM The Haag,
The Netherlands

Canada:

Health Canada device license number: 99490

Malaysian Medical Device Registration Number: GB10979720-50888

Malaysia Authorized Representative:
Emergo Malaysia Sdn. Bhd.
Level 16, 1 Sentral Jalan Stesen Sentral 5 KL Sentral, 50470
Kuala Lumpur MALAYSIA

CAUTION: *United States federal law restricts this device to sale by, or on the order of, a physician.*

Any serious incident that has occurred in relation to this device should be reported to NeoSoft and the competent authority of your Member State.

Safety

Introduction

To assure efficient and safe use it is essential to read this safety section and all associated topics before attempting to use the software. It is important for you to read and understand the contents of this manual before attempting to use this product. You should periodically review the procedures and safety precautions.

The software is intended for use by trained and qualified personnel only.

suiteDXT / suiteHEART® software has an expected useful service life of 7 years from its original release date.

CAUTION: Federal Law restricts this device to sale, distribution and use by or on the order of a physician.

The terms danger, warning, and caution are used throughout this manual to point out hazards and to designate a degree or level of seriousness. Hazard is defined as a source of potential injury to a person. Familiarize yourself with the terminology descriptions listed in the following table:

Table 1: Safety Terminology

Graphic	Definition
 DANGER:	Danger is used to identify conditions or actions for which a specific hazard is known to exist which <u>will</u> cause severe personal injury, death, or substantial property damage if the instructions are ignored.
 WARNING:	Warning is used to identify conditions or actions for which a specific hazard is known to exist.
 CAUTION:	Caution is used to identify conditions or actions for which a potential hazard is known to exist.

Indications for Use

suiteHEART® Software is an analytical software tool, which provides reproducible tools for the review and reporting of medical images. suiteHEART® Software can import medical images from a MR system and display them in a viewing area on the computer screen. The viewing area allows access to multiple studies and series of multi-slice, multi-phase images. Multi-phase sequences of images can be displayed in cine mode to facilitate visualization.

A report input interface is also available. Measurement tools on the report interface make it possible to quickly and reliably fill out a complete clinical report of an imaging exam. Available tools include: point, distance, area, and volume measurement tools such as ejection fraction, cardiac output, end-diastolic volume, end-systolic volume, and volume flow measurements.

Semi-automatic tools are available for left ventricular contour detection, valve plane detection, vessel contour detection for flow analysis, signal intensity analysis for myocardium and infarct sizing measurement, and T2 star analysis.

The results of the measurement tools are interpreted by the physician and can be communicated to referring physicians.

When interpreted by a trained physician these tools may be useful in supporting the determination of a diagnosis.

Intended Use

suiteHEART® Software is intended to assist trained clinical personnel in the qualification and quantification of cardiac function. The software provides the tools to adjust the parameters of the DICOM images and provides presentation states where the user can appreciate various MRI acquired images of the heart and vasculature over time. Additionally, the software provides tools for measuring linear distances, areas, and volumes that can be used to quantify cardiac function. Finally, the software provides the tools for volumetric flow measurements and the ability to calculate flow values.

Supported DICOM Image Formats

suiteHEART® Software supports the following DICOM format; MR and Enhanced MR. Refer to the suiteHEART® Software DICOM Conformance Statement manual for further detail on supported formats.

CAUTION: Data stored as a DICOM image that has been imported by an external PACS may not be compatible viewing for suiteHEART® Software.

Safety Notices

WARNING: The application assists in the analysis of the images only and does not automatically produce quantifiable results. The use and placement of quantitative measurements is at the discretion of the user. Misdiagnosis can occur if measurements are inaccurate. Measurements should only be created by a properly trained and qualified user.

WARNING: Artifacts on an image can be misinterpreted, leading to misdiagnosis. Do not use images containing artifacts for diagnosis. Analysis should only be accomplished by a properly trained and qualified user.

WARNING: Diagnosis for the wrong patient could occur if images do not contain patient name or ID. Do not use images that do not contain patient name and ID for diagnosis. Visually confirm the patient information prior to analysis.

CAUTION: Using images upon which an image filter has been applied could result in altered results. The user must exercise discretion before analyzing pixel intensity corrected images. The software will display a warning message if loading images that have been filtered.

Equipment Hazards

CAUTION: Using equipment that is damaged or has been compromised can put the patient at risk by delaying diagnosis. Make sure that equipment is in proper working order.

CAUTION: Applications may run on equipment that includes one or more hard disk drives, which may hold medical data related to patients. In some countries, such equipment may be subject to regulations concerning the processing of personal data and free circulation of such data. Release of personal data may result in legal action depending on the applicable regulatory body. It is strongly recommended that access to patient files be protected. The user is responsible for understanding the laws regulating patient information.

suiteDXT Installation & Configuration

Prerequisites

Prior to installation of the software, the following prerequisites must be met:

Table 1: System Requirements

Minimum System Requirements
Supported Operating Systems: • Windows 7 Professional or Enterprise with SP1 (64 bit version) • Windows 10 Professional or Enterprise (64 bit version)
Video card (NVIDIA); 4GB RAM, minimum of 1664 CUDA Cores, support for OpenGL 4.0, support 1920 x 1080 resolution or higher
Monitor shall support 1920 x 1080 resolution or higher
Anti-Virus/Anti-Malware software installed
Minimum 16GB memory
Available hard drive space minimum 50GB
PDF viewer, Adobe Reader 11.0 or higher
An open network port for DICOM transfer over network
Static IP address (recommended)
Solid-State hard drive
Intel® Core™ i7 processor
27" Monitor
Microsoft® .NET Framework 4.0

NOTE: In the event that the host system loses or is reassigned its IP address, remote devices will need to be manually reconfigured to request and accept network traffic.

NOTE: You must be logged into your operating system as Administrator in order to install the software, restart the DICOM service, and move the Local DICOM Storage.

Installation

1. Double click suiteDXT-windows-installer.
2. If prompted, select "Yes" on the User Account Control pop-up.
3. Select the desired language, then click "OK."

NOTE: If another version of the software is already installed on the system, a pop-up will ask permission to continue. Click "Yes" to continue with the installation.

4. Select "Next" on the suiteDXT setup screen.
5. Check the box for "I accept the agreement" and select "Next."
6. Uncheck the boxes for the software that is not being installed and select "Next."
7. Accept or change the installation directory and select "Next."
8. Select "Next" in ready to install setup screen.
9. Select "Finish" in the setup screen.

NOTE: DECISION POINT! Customer and NeoSoft should discuss enabling user level permissions.

Configuration

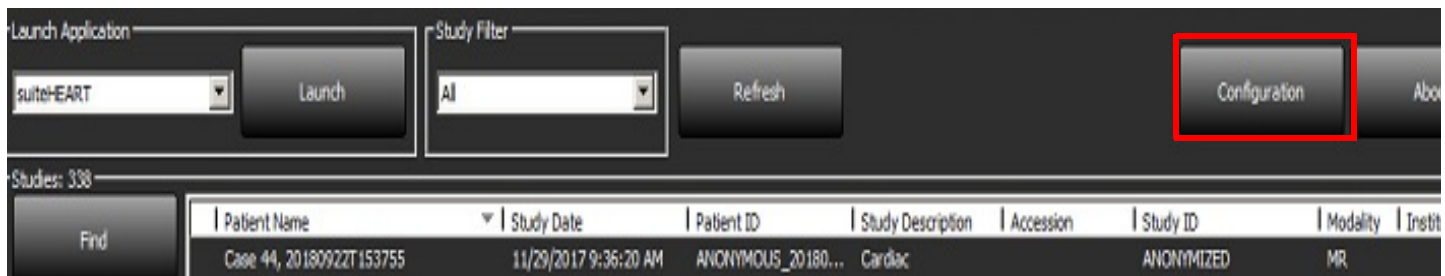
You will need to run through the configuration steps when suiteDXT is newly installed or when/if the static IP address is renewed.

1. Launch suiteDXT via the desktop shortcut or navigate through the start menu.

If prompted, select "Yes" on the user account pop-up.

2. Select the Configuration button, outlined in the figure below.

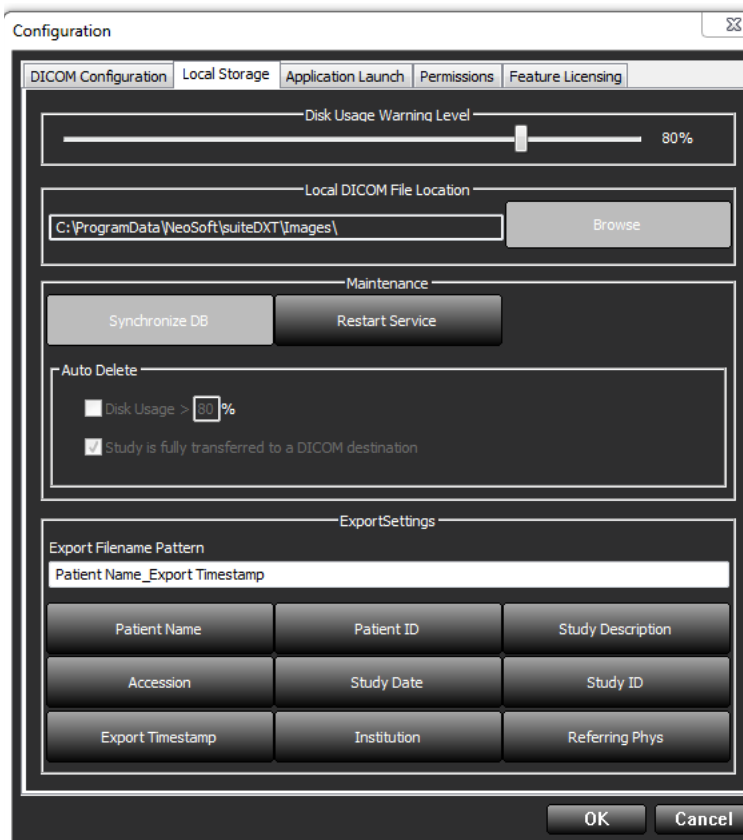
FIGURE 1. suiteDXT Main Screen

3. On the DICOM Configuration tab, enter information for the local host:
 - AE Title (case sensitive) is the name used to identify the host by other hosts.
 - Port (104).

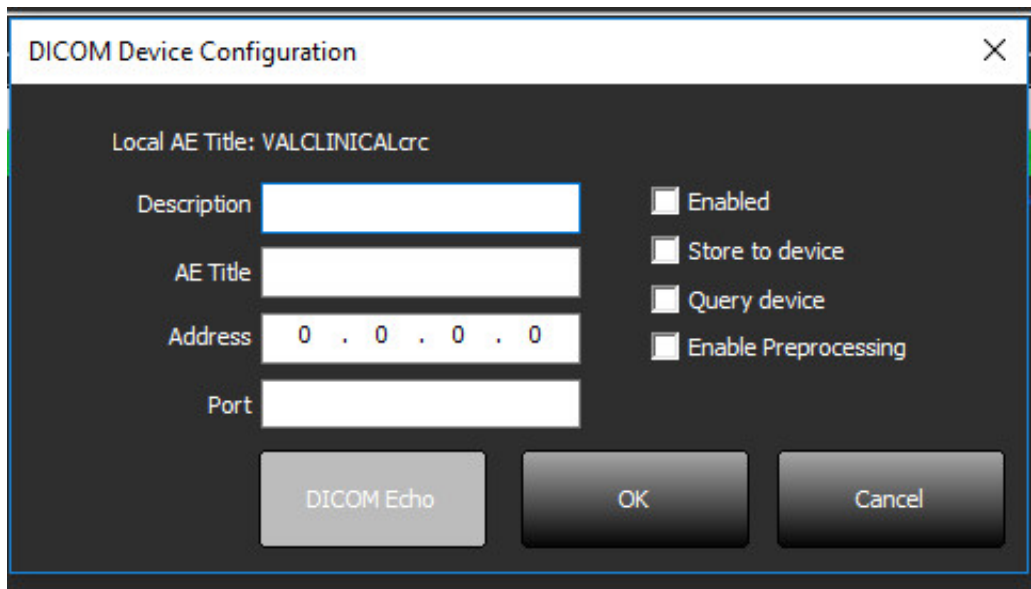
NOTE: The information entered on this tab must be entered exactly the same on any remote device for suiteDXT to retrieve or accept images from the remote host.

FIGURE 2. Configuration Tab Window

4. If available, perform a DICOM echo from any remote device to this host to ensure DICOM connectivity.
5. Select the “Add” button to input all DICOM destinations (remote hosts).
 - Description
 - AE Title
 - Address
 - Port
 - Select/deselect the checkboxes as necessary for: Enabled; Store to device, Query device and Enable Preprocessing.(Figure 3).

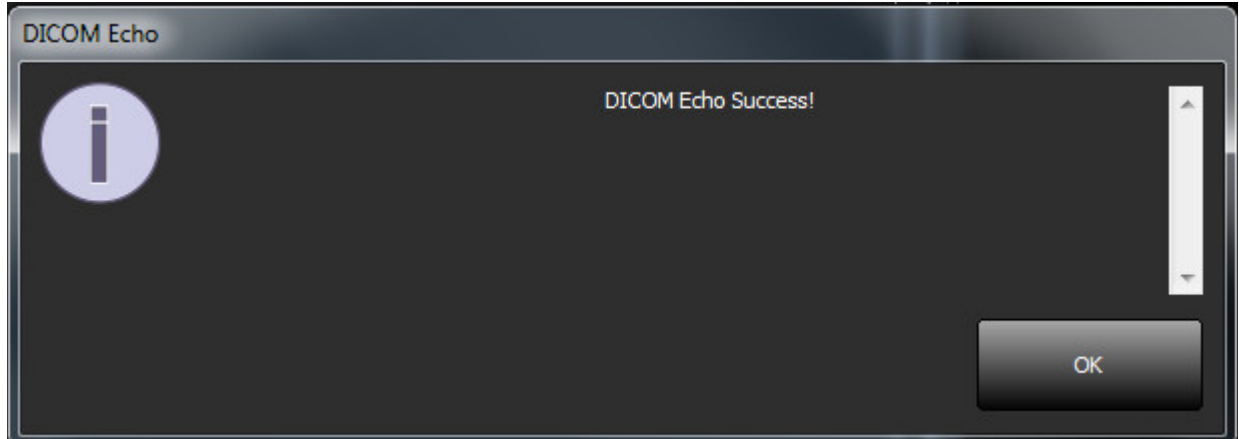
FIGURE 3. DICOM Configuration Window

The image shows a 'DICOM Device Configuration' window. At the top, it says 'Local AE Title: VALCLINICALcrc'. Below this are four input fields: 'Description', 'AE Title', 'Address' (with a default value of '0 . 0 . 0 . 0'), and 'Port'. To the right of these fields are four checkboxes: 'Enabled', 'Store to device', 'Query device', and 'Enable Preprocessing'. At the bottom of the window are three buttons: 'DICOM Echo', 'OK', and 'Cancel'.

- 6. Select the “DICOM Echo” button.
 - If you don’t see the DICOM Success Message window (Figure 4) contact your network administrator.

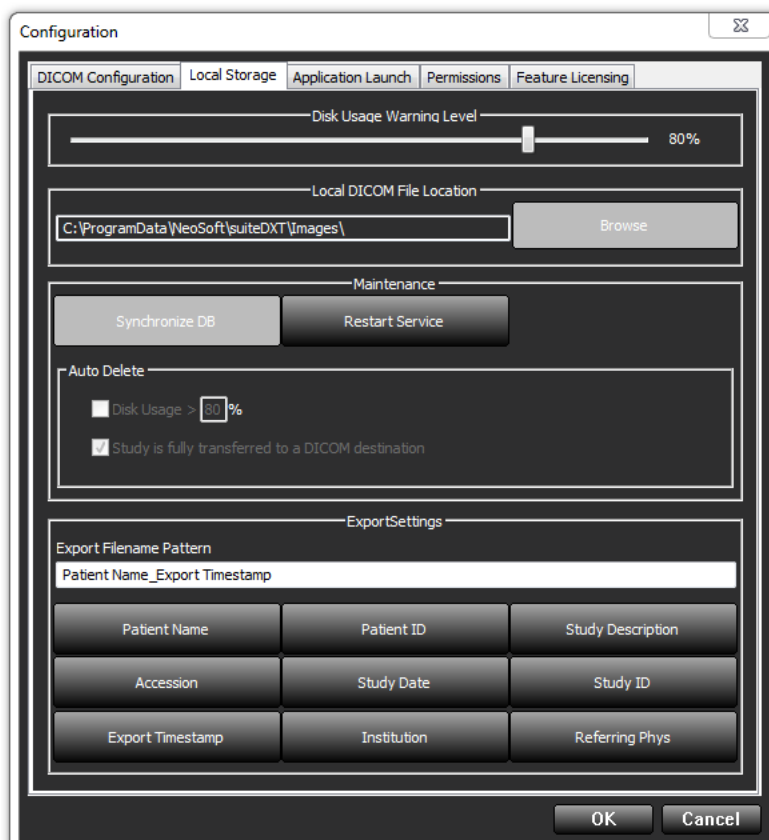
FIGURE 4. DICOM Success Message

NOTE: DICOM Destinations can be modified and a DICOM echo performed via the Edit button and destinations can be deleted via the Delete button.

Local Storage Tab

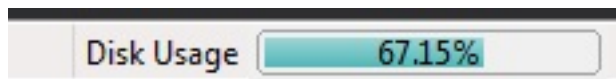
FIGURE 5. Local Storage Tab Window

The Disk Usage Warning Level can be modified to alert you when disk space is low. When disk usage exceeds the limit set, the indicator on the main screen (Figure 6) will change from blue to yellow.

To restore the study listing, click the Synchronize DB button. This refreshes the data currently stored on the local DICOM file location.

FIGURE 6. Disk Usage Indicator on Main Screen

The Local DICOM File Location Selection allows you to change the directory where DICOM data is stored. If data already exists it will be deleted before assigning a new folder. A warning message will appear to alert you to potential data loss.

CAUTION: Deleting DICOM data from the file system cannot be retrieved.

Auto Delete

The Auto Delete configuration should be maintained by the system administrator for your facility. Auto Delete can be configured to automatically delete the oldest studies in the system that match the defined criteria.

CAUTION: Studies that have been removed using the Auto Delete function cannot be retrieved.

Application Launch Tab

Use this tab to configure compatible applications to launch from suiteDXT.

1. Select the Add button (Figure 7) to open the Edit Application Window (Figure 8).

FIGURE 7. Application Launch Tab Window

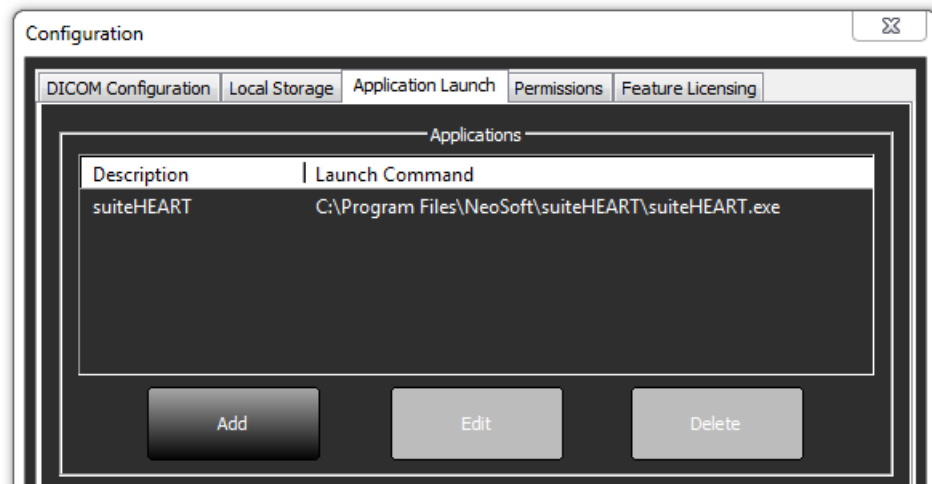
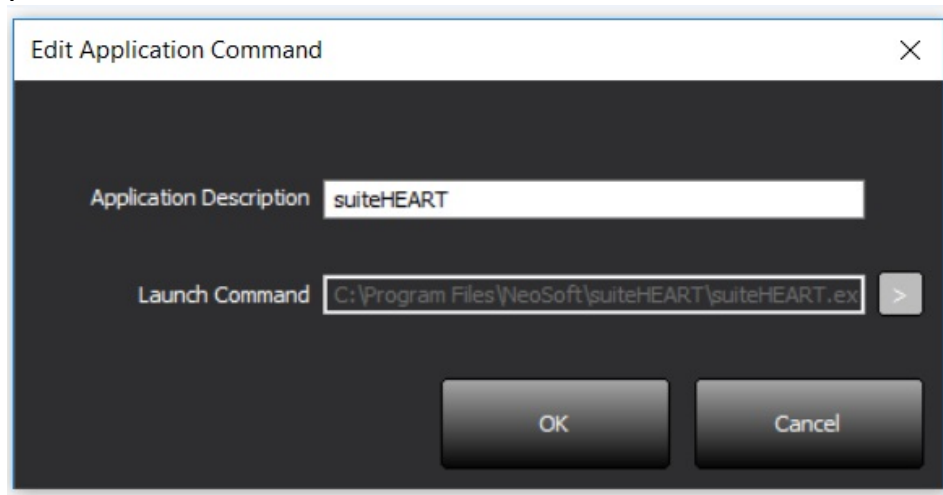
FIGURE 8. Edit Application Window

2. Enter the name or a description of the application.
3. Select the arrow button next to the Launch Command field to select the directory which contains the application executable file.
4. Select the OK button.

To edit or delete an existing application from the launch list, use the Edit or Delete button (Figure 7).

NeoSoft applications are automatically added into the application launch list. The launch command is not editable and nor can an application be manually deleted from the list.

Feature Licensing Tab

Use this tab to request and apply the license for all software distributed by NeoSoft, LLC.

FIGURE 9. Feature Licensing Tab Window

The screenshot shows a 'Configuration' window with a tabbed interface. The 'Feature Licensing' tab is selected. The window is divided into two main sections: 'Request License' and 'Apply License File'. The 'Request License' section contains two sub-sections: 'Company Information' with fields for Name, Telephone, and Billing Address; and 'Contact Information' with fields for First Name, Last Name, and Email. Below these fields is a 'Generate License Request' button. The 'Apply License File' section contains an 'Apply License' button. At the bottom right of the window are 'OK' and 'Cancel' buttons.

NOTE: Access to email will be needed to request a license.

Request a License

1. Enter complete company and contact information.

NOTE: Selecting the Generate License Request button will create a machine specific file that locks the software for which you are requesting a license, to the hardware (host) on which the License Server is installed.

2. Select the Generate License Request button.
3. Navigate to a memorable directory (e.g. the desktop) in which to save the license request file.
4. Email the license request file to licensing@neosoftmedical.com

Apply the License

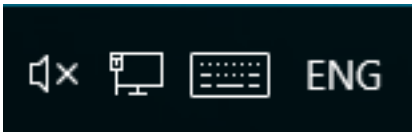
1. From NeoSoft's response email, copy the license to an easily accessible location (e.g. the desktop).
2. Open suiteDXT, select the Configuration button and navigate to the Feature Licensing tab (Figure 9).
3. Select the Apply License button and navigate to the directory to which you saved the license.
4. Select OK.

Tablet Mode

suiteDXT is supported on tablets with a 64-bit processor running on Windows 10 Professional or Windows 10 Enterprise operating system.

To activate the keyboard click on the windows keyboard icon located in the lower right corner.

FIGURE 10. Keyboard

suiteHEART® Software Installation

Prerequisites

Prior to installation of the software, the following prerequisites must be met:

Table 1: System Requirements

Minimum System Requirements
Supported Operating Systems: • Windows 7 Professional or Enterprise with SP1 (64 bit version) • Windows 10 Professional or Enterprise (64 bit version)
Video card (NVIDIA); 4GB RAM, minimum of 1664 CUDA Cores, support for OpenGL 4.0, support 1920 x 1080 resolution or higher
Monitor shall support 1920 x 1080 resolution or higher
Anti-Virus/Anti-Malware software installed
Minimum 16GB memory
Available hard drive space minimum 50GB
PDF viewer, Adobe Reader 11.0 or higher
An open network port for DICOM transfer over network
Static IP address (recommended)
Amazon Corretto version 15
Solid-State hard drive
Intel® Core™ i7 processor
27" Monitor

NOTE: suiteHEART® Software is not supported in Virtual Machine environments.

NOTE: suiteHEART® Software does not support scaling/zoom levels greater than 125%.

NOTE: You must be logged into your operating system as Administrator in order to install the software.

NOTE: For support of application frame rates of 30FPS at display resolutions in the range of 2K/4K, the NVIDIA video card hardware must support a minimum of 1920 CUDA Cores.

Installation

Software

1. Double-click on the suiteHEART® Software windows application installer.

NOTE: Installing the software on a machine could result in a Open File - Security Warning indicating that the publisher could not be verified. Click Yes when asked if you are sure you want to run the software.

2. Select the appropriate language for your installation and click **OK**.

NOTE: If another version of the software is already installed on the system, a pop-up will ask permission to uninstall. Click “Yes” to continue with the installation.

3. In the “Welcome to the suiteHEART® Setup Wizard” window, click **Next**.
4. On the License Agreement Window, select “I accept the agreement” and click **Next**.
5. Browse to the directory in which to install the software or use the default directory (recommended). Click **Next**.
6. Click **Next** on the Ready to Install screen.
A progress bar will display.
7. Click **Finish** to complete the installation.

Documentation

Prior to installing the documentation, the suiteDXT and suiteHEART® applications must be installed.

1. Double-click on the Documentation installer.

NOTE: Installing NeoSoft documentation on a machine could result in an Open File - Security warning indicating that the publisher could not be verified. Click **Yes** when asked if you are sure you want to run the software.

2. In the “Welcome to the NeoSoft Documentation Setup Wizard” window, click **Next**.
3. In the “Select Components” window, select the documentation to install (checking all available options is recommended), and click **Next**.
4. Click **Next** on the Ready to Install screen.
A progress bar will display.
5. Click **Finish** to complete the installation.

Anti-Virus Exclusions

To increase the performance of the application, the following Anti-Virus and Windows Defender exclusions should be made:

- C:\Program Files\NeoSoft, LLC or C:\Program Files\NeoSoft
- C:\ProgramData\NeoSoft, LLC or C:\ProgramData\NeoSoft

NOTE: The install locations noted above assume a default installation.

Application Inactivity Setting

The suiteHEART® Software application will automatically close, save all analysis for any open study and remove the associated lock file when the application has been inactive for 60 minutes. To set this setting configure it in suiteHEART® Software using Tools > Preferences > Edit.

Licensing Schemes

suiteHEART® Software has 3 different licensing schemes. They are Evaluation, Unlimited and Per-Case Pack Limited.

1. Evaluation: This is a time limited license which does not allow the saving of analysis results or exam approval. Evaluation Version will be indicated on the report footer and splash screens. If the approval of an exam is attempted a message will appear stating: Evaluation version: Analysis state will not be saved.
2. Unlimited: No restrictions on analysis and is not case-limited.
3. Per-Case Pack Limited: The software will use the unique study instance id to identify the number of cases analyzed against the Per-Case Pack that has been purchased.

Once 100% of the Per-Case Pack has been used the software will indicate this condition by a pop-message prior to starting a new exam or an existing exam.

suiteHEART® Software is licensed and sold using a user/seat license model.

Per-Case Packs can be purchased in increments of 25, 50 and 100. Website: <http://neosoftllc.com/> T: (262)522-6120
email: orders@neosoftmedical.com

Exams that have been counted once towards the Per-Case Pack limit can be re-opened and reviewed within the software and will not be counted against the Per-Case Pack limit.

Exams that are opened only for image viewing or if image manipulation tools are performed (Pan, Zoom, Window level, Rotate) will not be counted against the Per-Case Pack limit.

An exam is considered "consumed" or "counted" against the Per-Case Pack limit when any of the following actions are performed:

- a.) Starting of any analysis mode by placing any ROI on an image.
- b.) Creating a custom series.
- c.) Report sign off.
- d.) Export Cine DICOM.
- e.) Export Report.
- f.) DICOM series creation.
- g.) Preprocessed study.
- h.) Virtual Fellow® preprocessing.
- i.) Any Auto Combine / Split Mode.

suiteHEART® Software has an expected useful service life of 7 years from its original release date.

Remote Connectivity

Remote connectivity allows suiteHEART® Software users remote access from computers not loaded with suiteHEART® Software. The remote user will have the ability to access the suiteHEART® Software system and have full functionality.

suiteHEART® remote connectivity needs to be installed/configured/supported by your Information Technology personnel.

Remote connectivity has been verified using the following platforms:

- TeamViewer (version 15)
- Windows Remote Desktop (Microsoft Windows 7/10 Professional/Enterprise)

NOTE: When using Windows Remote Desktop, screen resolution on the remote computer must be set to 1920x1080 or higher video display resolution.

NOTE: When remotely accessing suiteHEART® Software using Windows Remote Desktop, the suiteHEART® 3D/4D viewer is not operational. The workaround is to use one of the other qualified remote access software packages (TeamViewer).

The performance of the suiteHEART® Software via the remote connection cannot be guaranteed. Performance is dependent upon factors that are outside the scope of the suiteHEART® Software application. These factors include:

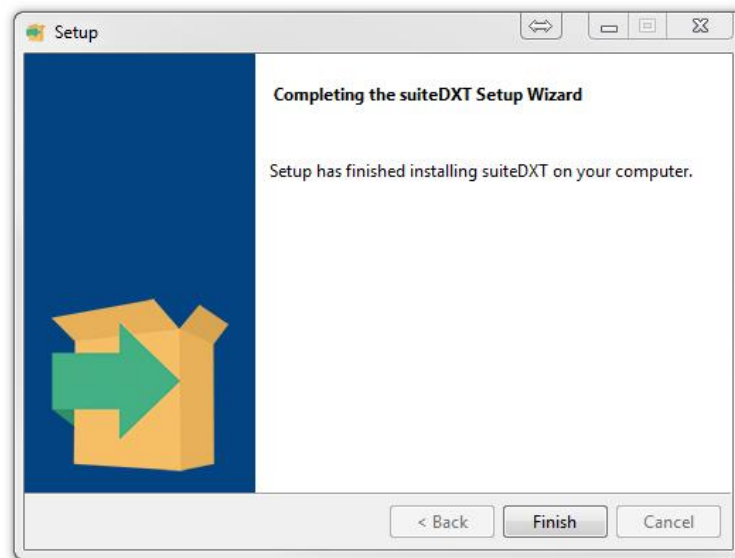
- Remote computer performance
- Internet download/upload speeds (Mbps)
- Network bandwidth speed (wired ethernet or wireless connection)

For more information on remote connectivity, refer to the suiteHEART® Software Addendum to Instructions for Use manual, NS-03-039-0014.

Install Verification

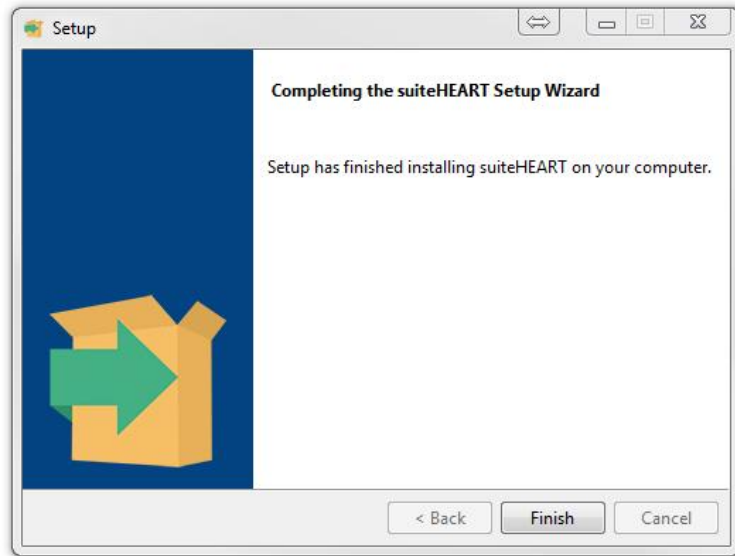
After a successful installation of suiteDXT this screen is displayed.

FIGURE 1. suiteDXT Setup Wizard

Following the successful installation of suiteHEART® Software this screen displays.

FIGURE 2. suiteHEART® Software SetupWizard

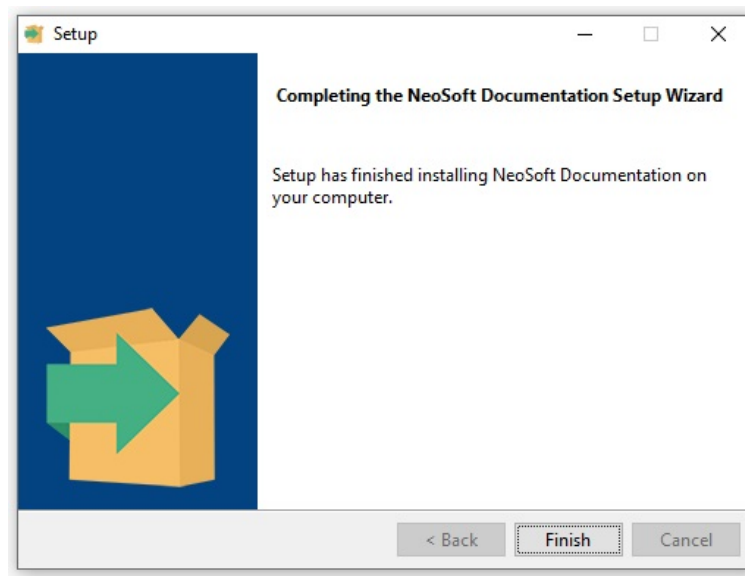
To further verify the successful installation of both software packages:

1. Double click the suiteDXT icon to verify suiteDXT is launched.
2. Highlight a study in the suiteDXT study list and select launch with suiteHEART to verify that the study is launched and that various analysis can be selected.

After a successful installation of the NeoSoft documentation this screen is displayed.

FIGURE 3. NeoSoft Documenatation Setup Wizard

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suiteHEART®-Software-Installation 12

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Malaysia Authorized Representative:

Emergo Malaysia Sdn. Bhd.
Level 16, 1 Sentral Jalan Stesen Sentral 5 KL Sentral, 50470
Kuala Lumpur MALAYSIA

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The image shows a software configuration window titled "Configuration". It has several tabs: "DICOM Configuration", "Local Storage", "Application Launch", "Permissions", and "Feature Licensing". The "Feature Licensing" tab is selected. Inside this tab, there is a "Request License" section. This section contains two sub-sections: "Company Information" and "Contact Information". "Company Information" has fields for "Name", "Telephone", and "Billing Address". "Contact Information" has fields for "First Name", "Last Name", and "Email". Below these fields is a button labeled "Generate License Request". Below the "Request License" section is an "Apply License File" section, which contains a button labeled "Apply License". At the bottom right of the window are "OK" and "Cancel" buttons.

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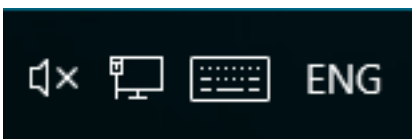
ABBILDUNG 10. Tastatur

suiteHEART®-Software-Installation

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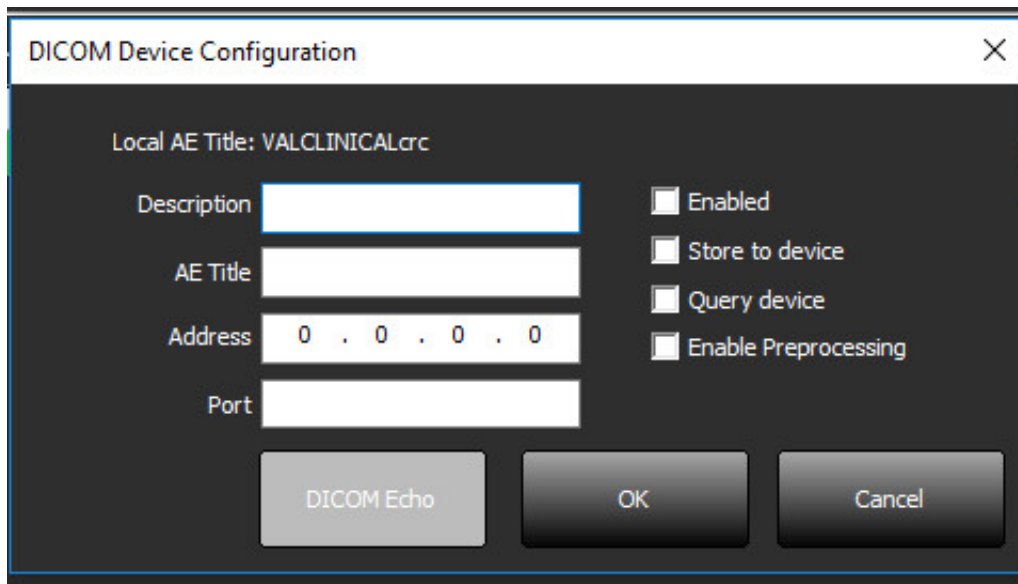
FIGURE 3. Fenêtre Configuration DICOM

The image shows a 'DICOM Device Configuration' dialog box. At the top, it says 'Local AE Title: VALCLINICALcrc'. Below this are four input fields: 'Description', 'AE Title', 'Address' (containing '0 . 0 . 0 . 0'), and 'Port'. To the right of these fields are four checkboxes: 'Enabled', 'Store to device', 'Query device', and 'Enable Preprocessing'. At the bottom of the dialog are three buttons: 'DICOM Echo', 'OK', and 'Cancel'.

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The image shows a software configuration window titled "Configuration". It has several tabs: "DICOM Configuration", "Local Storage", "Application Launch", "Permissions", and "Feature Licensing". The "Feature Licensing" tab is selected. Inside this tab, there are two main sections. The first section, titled "Request License", contains two sub-sections. The "Company Information" sub-section has input fields for "Name", "Telephone", and "Billing Address". The "Contact Information" sub-section has input fields for "First Name", "Last Name", and "Email". Below these fields is a button labeled "Generate License Request". The second section, titled "Apply License File", contains a button labeled "Apply License". At the bottom right of the window are "OK" and "Cancel" buttons.

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The image shows a software configuration window titled "Configuration". It has several tabs: "DICOM Configuration", "Local Storage", "Application Launch", "Permissions", and "Feature Licensing". The "Feature Licensing" tab is selected. Inside this tab, there is a section titled "Request License". This section contains two sub-sections: "Company Information" and "Contact Information". The "Company Information" sub-section has three text input fields: "Name", "Telephone", and "Billing Address". The "Contact Information" sub-section has three text input fields: "First Name", "Last Name", and "Email". Below these input fields is a button labeled "Generate License Request". Below this button is another section titled "Apply License File", which contains a button labeled "Apply License". At the bottom right of the window are "OK" and "Cancel" buttons.

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Configuration

DICOM Configuration | Local Storage | Application Launch | Permissions | Feature Licensing

Disk Usage Warning Level: 80%

Local DICOM File Location: C:\ProgramData\NeoSoft\suiteDXT\Images\ Browse

Maintenance: Synchronize DB Restart Service

Auto Delete: ☐ Disk Usage > 80% ☒ Study is fully transferred to a DICOM destination

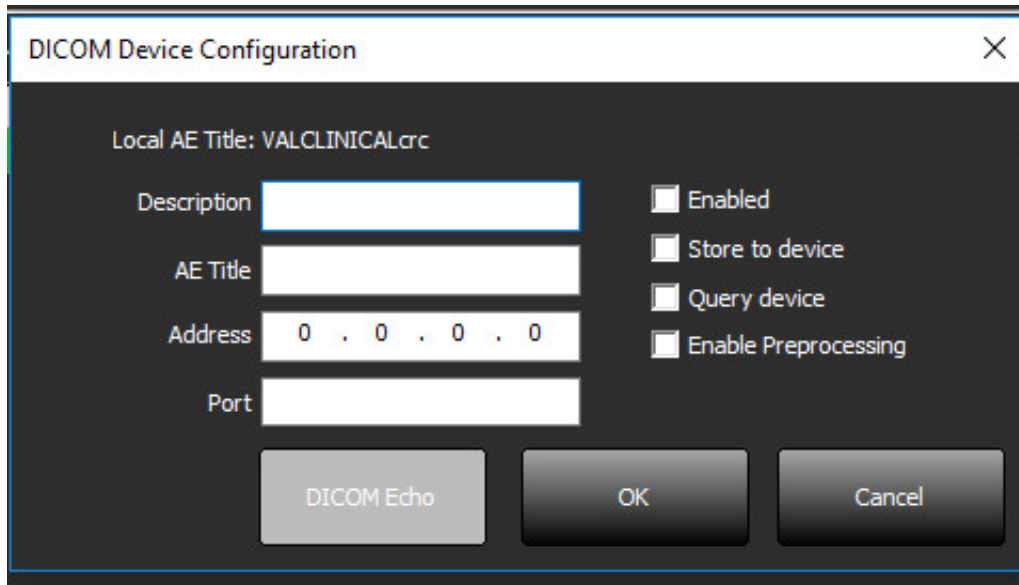
ExportSettings: Export Filename Pattern: Patient Name_Export Timestamp

Patient Name	Patient ID	Study Description
Accession	Study Date	Study ID
Export Timestamp	Institution	Referring Phys

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The image shows a 'DICOM Device Configuration' dialog box. At the top, it says 'Local AE Title: VALCLINICALcrc'. Below this are four input fields: 'Description', 'AE Title', 'Address' (with a default value of '0 . 0 . 0 . 0'), and 'Port'. To the right of these fields are four checkboxes: 'Enabled', 'Store to device', 'Query device', and 'Enable Preprocessing'. At the bottom of the dialog are three buttons: 'DICOM Echo', 'OK', and 'Cancel'.

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The image shows a software configuration window titled "Configuration". It has several tabs: "DICOM Configuration", "Local Storage", "Application Launch", "Permissions", and "Feature Licensing". The "Feature Licensing" tab is selected. Inside this tab, there is a section titled "Request License". This section contains two sub-sections: "Company Information" and "Contact Information". The "Company Information" sub-section has three text input fields: "Name", "Telephone", and "Billing Address". The "Contact Information" sub-section has three text input fields: "First Name", "Last Name", and "Email". Below these fields is a button labeled "Generate License Request". Below the "Request License" section is another section titled "Apply License File", which contains a button labeled "Apply License". At the bottom right of the window are two buttons: "OK" and "Cancel".

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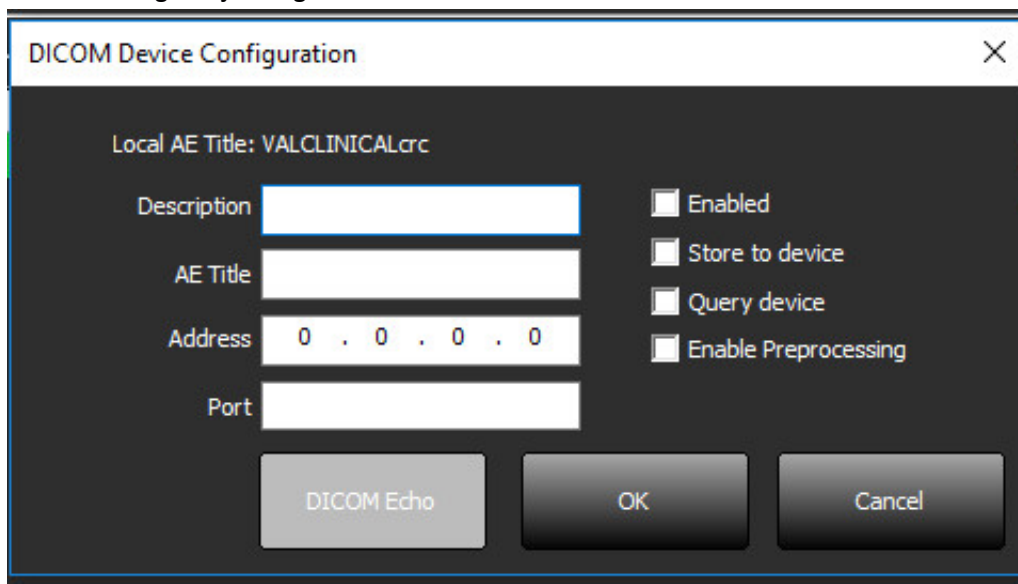
The image shows a 'DICOM Device Configuration' dialog box. At the top, it says 'Local AE Title: VALCLINICALcrc'. Below this, there are four input fields: 'Description', 'AE Title', 'Address' (containing '0 . 0 . 0 . 0'), and 'Port'. To the right of these fields are four checkboxes: 'Enabled', 'Store to device', 'Query device', and 'Enable Preprocessing'. At the bottom, there are three buttons: 'DICOM Echo', 'OK', and 'Cancel'.

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The screenshot shows a 'Configuration' window with the 'Feature Licensing' tab selected. The window contains two main sections: 'Request License' and 'Apply License File'. The 'Request License' section includes 'Company Information' (Name, Telephone, Billing Address) and 'Contact Information' (First Name, Last Name, Email). Below these is a 'Generate License Request' button. The 'Apply License File' section contains an 'Apply License' button. At the bottom of the window are 'OK' and 'Cancel' buttons.

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The screenshot shows a 'Configuration' window with a tabbed interface. The 'Feature Licensing' tab is selected. The window is divided into two main sections: 'Request License' and 'Apply License File'. The 'Request License' section contains two sub-sections: 'Company Information' with fields for Name, Telephone, and Billing Address; and 'Contact Information' with fields for First Name, Last Name, and Email. Below these is a 'Generate License Request' button. The 'Apply License File' section contains an 'Apply License' button. At the bottom right are 'OK' and 'Cancel' buttons.

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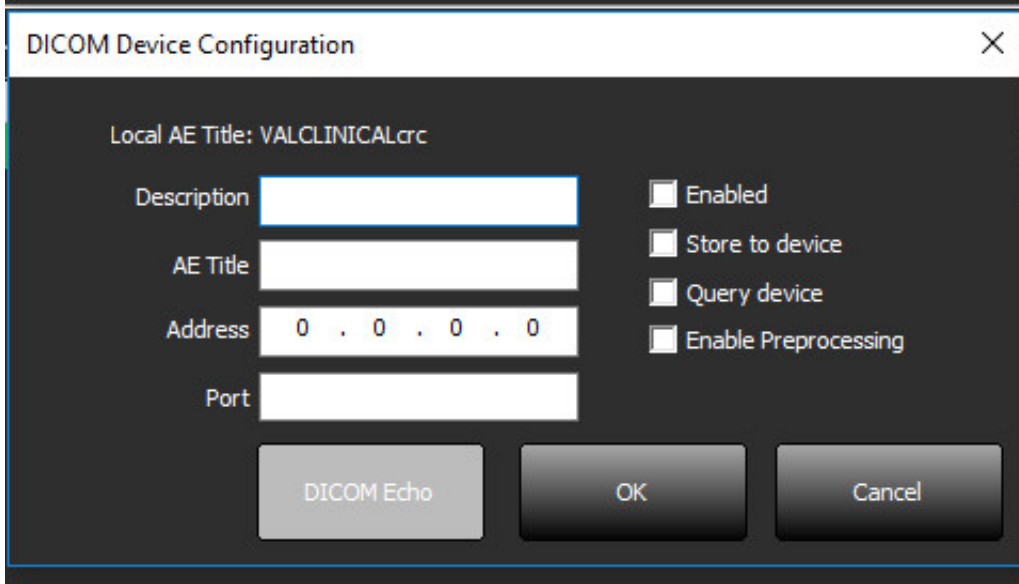
The Configuration window has a title bar with a close button. It contains several tabs: DICOM Configuration, Local Storage, Application Launch, Permissions, and Feature Licensing. The DICOM Configuration tab is active. It features a slider for Disk Usage Warning Level set at 80%. Below this is a text field for Local DICOM File Location with the path C:\ProgramData\NeoSoft\suiteDXT\Images\ and a Browse button. The Maintenance section has buttons for Synchronize DB and Restart Service. The Auto Delete section has checkboxes for Disk Usage > 80% (unchecked) and Study is fully transferred to a DICOM destination (checked). The ExportSettings section has a text field for Export Filename Pattern with the value Patient Name_Export Timestamp, and a grid of buttons for Patient Name, Patient ID, Study Description, Accession, Study Date, Study ID, Export Timestamp, Institution, and Referring Phys. At the bottom are OK and Cancel buttons.

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The image shows a 'DICOM Device Configuration' window. At the top, it says 'Local AE Title: VALCLINICALcrc'. Below this, there are four input fields: 'Description', 'AE Title', 'Address' (with a default value of '0 . 0 . 0 . 0'), and 'Port'. To the right of these fields are four checkboxes: 'Enabled', 'Store to device', 'Query device', and 'Enable Preprocessing'. At the bottom of the window are three buttons: 'DICOM Echo', 'OK', and 'Cancel'.

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The screenshot shows a 'Configuration' window with the 'Feature Licensing' tab selected. The window is divided into two main sections: 'Request License' and 'Apply License File'. The 'Request License' section contains two sub-sections: 'Company Information' with fields for Name, Telephone, and Billing Address; and 'Contact Information' with fields for First Name, Last Name, and Email. Below these fields is a 'Generate License Request' button. The 'Apply License File' section contains an 'Apply License' button. At the bottom of the window are 'OK' and 'Cancel' buttons.

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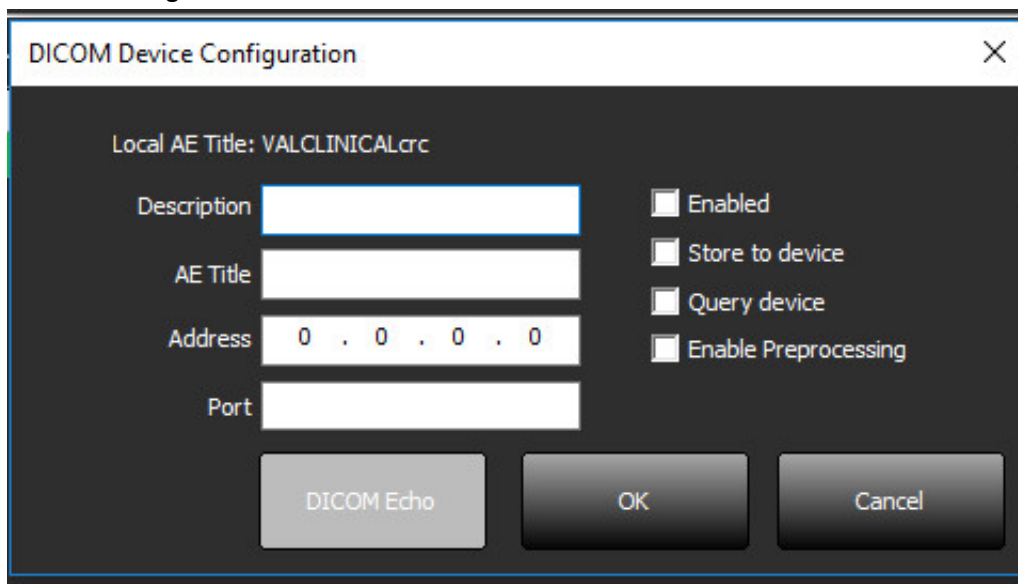
FIGURA 3. Fereastra de configurare DICOM

The image shows a 'DICOM Device Configuration' dialog box. At the top, it says 'Local AE Title: VALCLINICALcrc'. Below this are four input fields: 'Description', 'AE Title', 'Address' (with a dotted pattern 0 . 0 . 0 . 0), and 'Port'. To the right of these fields are four checkboxes: 'Enabled', 'Store to device', 'Query device', and 'Enable Preprocessing'. At the bottom, there are three buttons: 'DICOM Echo', 'OK', and 'Cancel'.

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The image shows a 'Configuration' window with the 'Feature Licensing' tab selected. The window contains two main sections: 'Request License' and 'Apply License File'. The 'Request License' section includes 'Company Information' (Name, Telephone, Billing Address) and 'Contact Information' (First Name, Last Name, Email). Below these is a 'Generate License Request' button. The 'Apply License File' section contains an 'Apply License' button. At the bottom of the window are 'OK' and 'Cancel' buttons.

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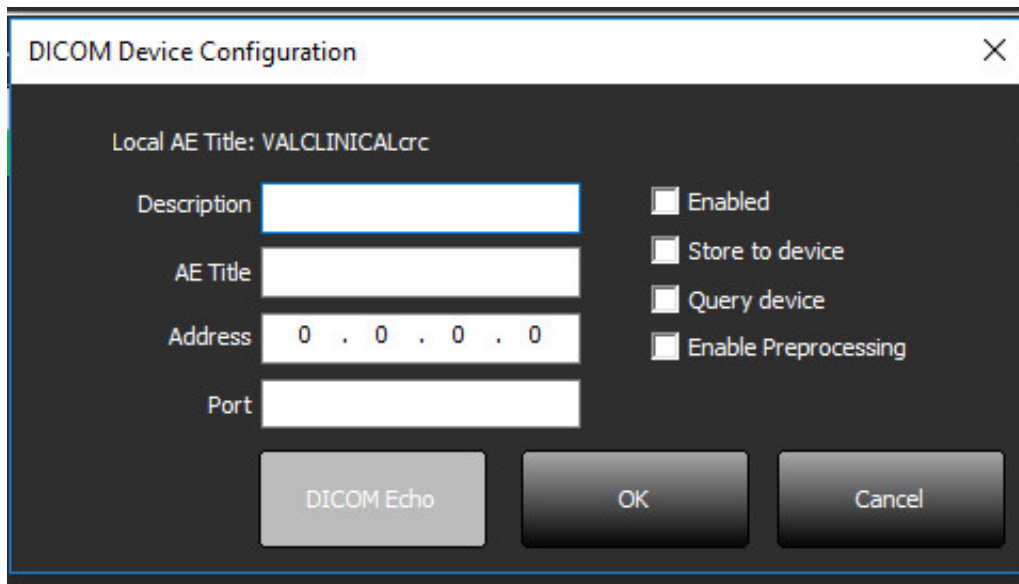
The image shows a 'DICOM Device Configuration' window. At the top, it says 'Local AE Title: VALCLINICALcrc'. Below this are four input fields: 'Description', 'AE Title', 'Address' (containing '0 . 0 . 0 . 0'), and 'Port'. To the right of these fields are four checkboxes: 'Enabled', 'Store to device', 'Query device', and 'Enable Preprocessing'. At the bottom of the window are three buttons: 'DICOM Echo', 'OK', and 'Cancel'.

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The screenshot shows a 'Configuration' window with a tabbed interface. The 'Feature Licensing' tab is selected. The window is divided into two main sections: 'Request License' and 'Apply License File'. The 'Request License' section contains two sub-sections: 'Company Information' with fields for 'Name', 'Telephone', and 'Billing Address'; and 'Contact Information' with fields for 'First Name', 'Last Name', and 'Email'. Below these fields is a 'Generate License Request' button. The 'Apply License File' section contains an 'Apply License' button. At the bottom right of the window are 'OK' and 'Cancel' buttons.

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FIGURA 9. Janela de separador de Licenciamento de recursos

The image shows a 'Configuration' window with a tabbed interface. The 'Feature Licensing' tab is selected. Inside this tab, there are two main sections: 'Request License' and 'Apply License File'. The 'Request License' section contains two sub-sections: 'Company Information' with fields for 'Name', 'Telephone', and 'Billing Address'; and 'Contact Information' with fields for 'First Name', 'Last Name', and 'Email'. Below these fields is a 'Generate License Request' button. The 'Apply License File' section contains an 'Apply License' button. At the bottom right of the window are 'OK' and 'Cancel' buttons.

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The image shows a software configuration window titled "Configuration". It has several tabs: "DICOM Configuration", "Local Storage", "Application Launch", "Permissions", and "Feature Licensing". The "Feature Licensing" tab is selected. Inside this tab, there is a "Request License" section. This section contains two sub-sections: "Company Information" and "Contact Information". "Company Information" has fields for "Name", "Telephone", and "Billing Address". "Contact Information" has fields for "First Name", "Last Name", and "Email". Below these fields is a button labeled "Generate License Request". At the bottom of the "Request License" section is another button labeled "Apply License". At the very bottom of the dialog are "OK" and "Cancel" buttons.

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